

TRIARCH

ACCURACY

DEPENDABILITY

SERVICE

Are Your Pet Slides Fading?

Preserve their teaching value by having them photomicrographed. Lantern slides from these photomicrographs help wonderfully in class presentation.

Triarch photomicrographic prints are like fine pen-and-ink drawings—clean, white backgrounds—just right for snappy half-tone reproduction.

*Let Us Make Your Photomicrographs and
Lantern Slides: Triarch Experts Know How!*

TRIARCH BOTANICAL PRODUCTS


"Prepared by a Botanist for Botanists"

include:


PREPARED SLIDES FOR GENERAL BOTANY. SPECIAL PREPARATIONS FOR PLANT PATHOLOGY. PRESERVED PLANTS. ALGAE IN PURE CULTURE. BOTANICAL MODELS. RESEARCH SERVICE. MICROSCOPES. MICROPROJECTORS. LABORATORY EQUIPMENT. PHOTOMICROGRAPH AND LANTERN SLIDE SERVICE.

For catalogs and information write:

GEO. H. CONANT . . . **Ripon, Wisconsin**

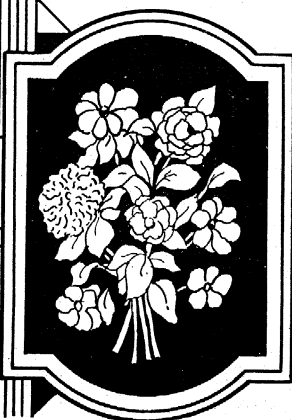

BOTANICAL PRODUCTS

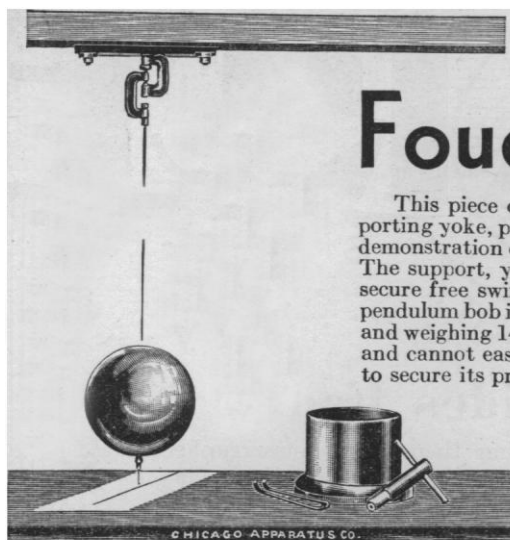

Foucault Pendulum

This piece of apparatus consisting of a massive pendulum bob, supporting yoke, pivot, wire, etc., is especially designed for the clear, effective demonstration of the Earth's rotation and for the study of the pendulum. The support, yokes and bearings, pivot, etc., are designed and fitted to secure free swinging action but safety and ease in demonstration. The pendulum bob is an accurately turned steel sphere about 14 cm. in diameter and weighing 14.5 kilos. It is so constructed as to secure homogeneous mass and cannot easily be distorted by accidents. It carries fittings necessary to secure its proper performance in experiments.


This apparatus is furnished with facilities making it easy to handle the ball for storage and transport. Special wrenches are also supplied. A storage support for the ball is shown in the illustration.

Price complete - - - - - \$40.00



CHICAGO APPARATUS COMPANY

1735-43 No. Ashland Avenue, Chicago, Ill.

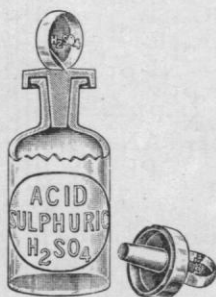

1894

Reduced Prices for JEWETT REAGENT BOTTLES

These bottles were designed by Professor F. F. Jewett, of Oberlin College, Oberlin, Ohio, and have been well received by those desiring a really practical Reagent Bottle having not only a protection against dirt, but a protection against the cementation of stoppers in the bottles, especially those containing alkalis.

The number of these bottles now being purchased enables us to materially reduce the price, and we hope that this inducement will stimulate the sale of these bottles generally.

FOUR OUNCE SIZE		REDUCED PRICES		EIGHT OUNCE SIZE	
Each		\$0.60		Each	\$0.70
Per Dozen		6.00		Per Dozen	7.00

These bottles may be obtained with raised letters in a wide assortment of names. If, however, it is desired to have bottles with other names than those contained in our printed lists, names with letters to match may be engraved on the bottles. The difference in the appearance of the raised letters and engraved letters is not sufficient to attract attention.


Printed lists or the available names as carried in stock cheerfully furnished upon request



E.H. SARGENT & CO.


155-165 East Superior Street · · · CHICAGO, ILL.
LABORATORY SUPPLIES

(4092B)

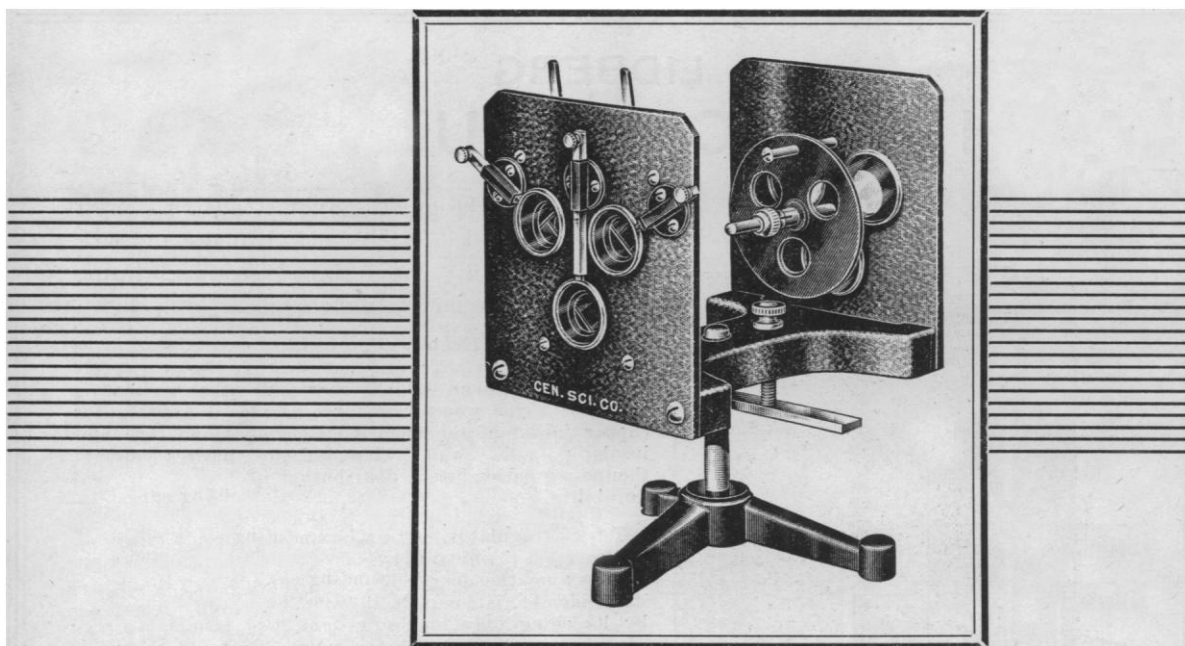


For the Study of Color Phenomena «

The Von Nardroff Color Apparatus «

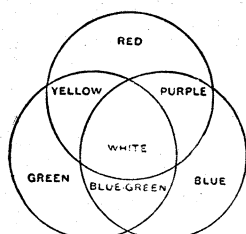

MAXWELL'S PRIMARY HUES.
Fig. 1.

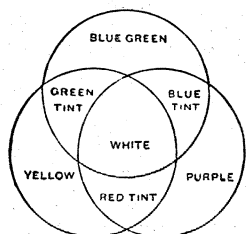

FALSE PRIMARY HUES
Fig. 2.


The Von Nardroff Color Apparatus was produced for the convenient and satisfactory study of color phenomena through the use of only one projecting lantern or other intense source of illumination. The single beam of light from the illuminator is divided into three independent beams of light, each of which can be modified in intensity or color by means of diaphragms and color screens. By moving the lens arms these modified beams of light may be overlapped either in whole or in part as desired and by adjusting the diaphragms the intensity or strength of each color can be regulated to produce a wide variety of shades in the mixed colors. Shadow effects can be obtained by holding an opaque object in the path of the light beams.

The Von Nardroff Color Apparatus is furnished complete with six color filters and directions for use.

No. F7175 \$38.00

CENTRAL SCIENTIFIC COMPANY
LABORATORY SUPPLIES
Apparatus and Chemicals
NEW YORK - BOSTON - CHICAGO - TORONTO - LOS ANGELES

LIDBERG ANHYDRIC INCUBATORS

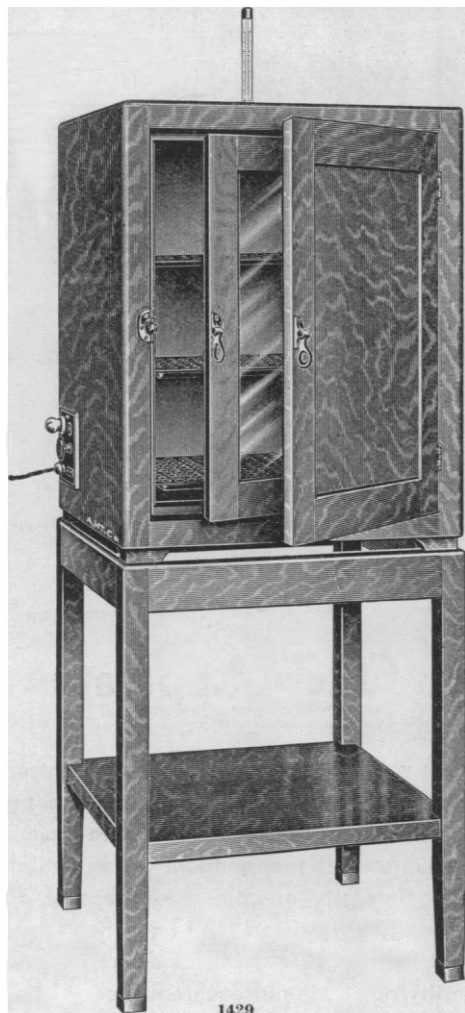

A newly designed anhydric construction of 5-ply laminated quartered oak, transite and hair felt insulation, with tinned copper lining throughout the incubating chamber.

Woven Mesh Heater. The heating unit consists of a woven resistance wire and asbestos mesh—an entirely novel arrangement—which surrounds the incubating chamber. This woven resistance is placed *between* the copper jacket of the incubating chamber and the outer insulated wall, which arrangement insures prompt thermo-regulation, even distribution of heat, and great durability.

Thermo-regulator. Of an expanding capsule type, controlled from the outside by a fine adjustment screw without opening the door. The regulating mechanism, with its condenser, is complete in one unit, readily removable by loosening only one screw and two terminals. The capsule operates through the range 36° to 46°C, and is quickly replaceable. It operates equally well on a.c. or d.c. and, since it is entirely outside the incubating chamber, no sparking contacts occur within the chamber.

Shelf Supports. Adjustable for the height best suited for special purposes or standardized loads. The shelves are of perforated metal and slide in runners for ready withdrawal of contents.

Entire capacity of incubating chamber available. No space is taken up by heaters, thermo-regulators, or unusable non-isothermic areas.

Operating sensitivity. Within $\pm 0.5^\circ\text{C}$ at 37°C.

Operating uniformity. Within $\pm 2.0^\circ\text{C}$ at 37°C.

Unit system. Made in 18-unit and 48-unit sizes.

Current consumption. No. 1423 consumes 140 watts, and No. 1428 consumes 220 watts on either a.c. or d.c. lighting circuit of 110 volts.

Floor Space. The 18-unit size requires a floor or table space of 21½ × 21½ inches, and the 48-unit size a space of 26 × 26 inches.

		Code Word
1423.	Incubator, Lidberg Anhydric, 18-unit, as above described, with net available inside dimensions 15½ inches wide × 15½ inches deep × 18 inches high (4,324 cubic inches). With thermometer with 180 mm stem and complete operating directions. On low stand, i.e. wooden feet 1-inch high, for use on table. For 110 volts	180.00 Aibai
1427.	Ditto, but on high stand, as shown in illustration	200.00 Aibea
1428.	Incubator, Lidberg Anhydric, 48-unit, as above described, with net available inside dimensions 20 inches wide × 20 inches deep × 26½ inches high (10,600 cubic inches). With thermometer with 180 mm stem and complete operating directions. On low stand, i.e. wooden feet 1-inch high, for use on table. For 110 volts	225.00 Aibko
1429.	Ditto, but on high stand, as shown in illustration	250.00 Aibmk

Prices and specifications subject to change without notice.

ARTHUR H. THOMAS COMPANY

RETAIL—WHOLESALE—EXPORT

LABORATORY APPARATUS AND REAGENTS

WEST WASHINGTON SQUARE

Cable Address, BALANCE, Philadelphia

PHILADELPHIA, U. S. A.

SCIENCE

VOL. 74

FRIDAY, AUGUST 7, 1931

No. 1910

<i>Observations upon the Filterability of Bacteria, including a Filterable Organism Obtained from Cases of Influenza:</i> PROFESSOR ARTHUR ISAAC KENDALL	129	<i>Societies and Academies:</i> <i>The North Carolina Academy of Science:</i> DR. H. R. TOTTEN. <i>The Virginia Academy of Science:</i> DR. E. C. L. MILLER	153
<i>Arrest of Geologic, Archeologic and Paleontologic Work in Central Asia:</i> PRESIDENT HENRY FAIRFIELD OSBORN	139	<i>Scientific Apparatus and Laboratory Methods:</i> <i>Ruthenium Tetroxide as a Fixative in Cytology:</i> DR. D. C. CARPENTER and DR. B. R. NEBEL. <i>Method for Locating the Larvae of the Mosquito <i>Mansonia</i>:</i> T. E. MCNEEL	154
<i>Obituary:</i> <i>Aldred Scott Warthin:</i> DR. PEYTON ROUS. <i>Memorials; Recent Deaths</i>	142	<i>Special Articles:</i> <i>On Atmospheric Electricity:</i> PROFESSOR C. S. DORCHESTER and PROFESSOR L. W. BUTLER. <i>Use of Protozoa in Measuring the Neutralizing Value of Cobra Anti-Serum:</i> DR. CHARLES H. PHILPOTT... ..	155
<i>Scientific Events:</i> <i>The Jubilee Meeting of the British Society of Chemical Industry; The Buckston-Browne Surgical Research Farm; The Sixth International Congress of Genetics</i>	144	<i>Science News</i>	10
<i>Scientific Notes and News</i>	146		
<i>Discussion:</i> <i>Reduction of Oxyluciferin by Atomic Hydrogen:</i> PROFESSOR E. NEWTON HARVEY and DR. G. I. LAVIN. <i>"The Possible Rôle of Micro-organisms in the Precipitation of Calcium Carbonate in Tropical Seas":</i> DR. CHAS. B. LIPMAN. <i>Inherited Taste Deficiency:</i> DR. L. H. SNYDER. <i>Getting the Student to Use His Own Intellect:</i> DR. M. LOUISE NICHOLS. <i>Positive Gas Pressure in Poplar:</i> PROFESSOR FRANK C. GATES	150		

SCIENCE: A Weekly Journal devoted to the Advancement of Science, edited by J. McKEEN CATTELL and published every Friday by

THE SCIENCE PRESS

New York City: Grand Central Terminal
Lancaster, Pa. Garrison, N. Y.
Annual Subscription, \$6.00 Single Copies, 15 Cts.

SCIENCE is the official organ of the American Association for the Advancement of Science. Information regarding membership in the Association may be secured from the office of the permanent secretary, in the Smithsonian Institution Building, Washington, D. C.

OBSERVATIONS UPON THE FILTERABILITY OF BACTERIA, INCLUDING A FILTERABLE ORGANISM OBTAINED FROM CASES OF INFLUENZA¹

STUDIES IN BACTERIAL METABOLISM, CI


By Professor ARTHUR ISAAC KENDALL

DEPARTMENT OF RESEARCH BACTERIOLOGY, NORTHWESTERN UNIVERSITY MEDICAL SCHOOL,
CHICAGO, ILLINOIS

THERE is a heterogenous group of formidable diseases of man and animals, including some of the most rapidly spreading infections, whose inciting agents have thus far eluded cultivation upon artificial mediums. There can be little doubt that the incitants of this group are living entities and their cultivation outside the body should, judging from past experi-

ence in bacteriology and preventive medicine, be helpful, not only in solving important problems of their life history, but also in approaching the solution of prophylactic and curative measures. A few of these "viruses" are said to have been kept alive for periods of time in presence of large amounts of blood, or pieces of tissue from recently killed animals, but by common consent, this limited, restricted development is not regarded as equivalent to cultivation in the usual bacteriological sense.

¹ The James A. Patten Lecture in Bacteriology, Northwestern University Medical School, Chicago, Illinois, July 22, 1931.

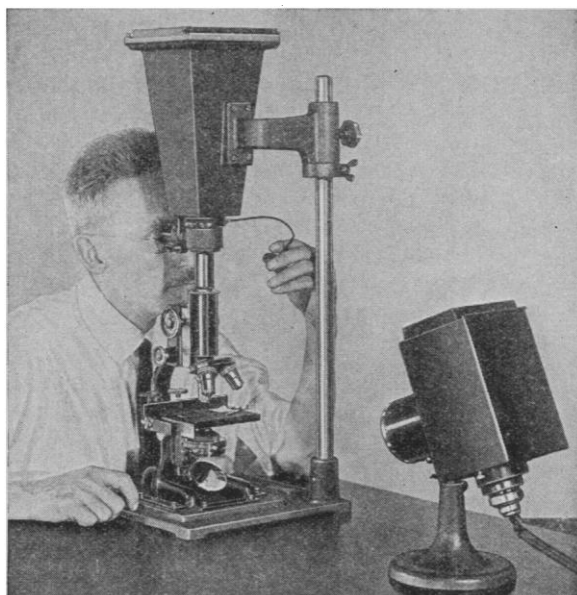

When You Make Photomicrographs

THE Model K Camera has been developed for routine photomicrography. It is simply constructed for use with any standard microscope and light source. An observation eyepiece permits the specimen to be examined and photographed without altering the position of either the microscope or the camera, and the specimen may be observed while the exposure is being made.

The camera body is supported entirely independently of the

microscope and its weight does not rest on the fine adjustment.

Model K is an all metal camera. It is adjustable upon a steel rod, which is mounted upon a cast iron base. The camera back is of standard form, and takes a $3\frac{1}{4} \times 4\frac{1}{4}$ double plate holder, a cut film holder, a film pack adapter or a roll film holder.

Write for complete information. BAUSCH & LOMB OPTICAL COMPANY, 642 St. Paul Street, Rochester, New York.



BAUSCH & LOMB

McGraw-Hill

has just published

PHYSICAL SCIENCE

An Introduction to the Specialized
Courses in College Science

By Herbert Brownell

*Professor of the Technique of Instruction
in Science, Teachers College, University of Nebraska*

313 pages, 6x9, 125 illustrations

\$2.50

The author has prepared this text with the following objectives in view:

- to present in a one-semester's course of three college hours credit enough of suitably related subject matter from the different fields of physical science to serve as an **Orientation Course** in science, for high school graduates entering upon college work, and as a summarization and reorganization of whatever science has been gained by them in high school.
- to make use of sufficient laboratory exercises, chiefly demonstration experimental work and field observations, to maintain the spirit and educational aims of science instruction, thus making it a **science course** as distinguished from a course **about science**.
- to ignore the dividing lines between the various fields of specialized science, but at the same time through choices of topics to maintain a suitable representation and balance of subject-matter from all the fields.
- to project into a post-high school stage of school work an attempt to interpret environment in ways suitably characterized, it may be, as **what everyone should know** of the teachings of science as these directly affect our daily lives.

Send for a copy on approval

McGRAW-HILL BOOK COMPANY, Inc.

Penn Terminal Building

370 Seventh Avenue

New York

*Books that have a Significant Place in the Library of
the Engineer, the Chemist and Technologist*

DIELECTRIC CONSTANT and MOLECULAR STRUCTURE

Charles Phelps Smyth

224 Pages \$4.00

A. C. S. Monograph No. 55

Just Published

This monograph collects, correlates, and interprets a mass of material which has recently developed in a new field of very general interest. The modern mathematical theory of the dielectric constant is outlined and its use in calculating the electric moments of molecules is shown. Phenomena related to dielectric constant and moment are described, among the subjects treated being molecular association and anomalous dispersion for radio frequencies, which is of special importance in problems of insulation. The main body of the book is devoted to the general relation of the electric or dipole moment to the structure of the molecule and to the particular relation of moment to structure in many molecules, radicals, and classes of compounds. Such problems as the structure of the benzene ring, electronic theories of valence, and the relation of refraction to structure are treated. The book makes available to the physicist, the chemist, and the electrical engineer the results of many recent investigations, showing the success of a new method of attack upon the molecule and describing methods for the further prosecution of the attack.

Chapter Headings: 1. Mathematical Theory of the Dielectric Constant. 2. Phenomena Related to Dielectric Constant and Electric Moment. 3. The Measurement of Dielectric Constant. 4. The Relation of Electric Moment to Molecular Structure. 5. The Electric Moments of Molecules and Radicals. 6. The Electric Moments and Molecular Structures of Cyclic Compounds. 7. Electric Moments and Electronic Theories of Valence. 8. Electronic and Atomic Polarization. 9. Molecular Association. Appendix. Table of Electric Moments with References. With Thirty-one Diagrams.

This work brings together the various procedures and applications of the volumetric iodate methods employed in routine analysis and research work.

Contents: The iodate method; The titration of iodine, iodides, etc.; The preparation of standard potassium iodate solutions; The determination of antimony; Antimony in alloys; Arsenic; Arsenic in insecticides; Copper; Copper in ores; Copper in alloys; Copper in insecticides; Mercury (mercuric); Mercury in antiseptic tablets; mercurous chloride; The titration of mercurous chloride in tablets; The determination of molybdenum; Tin; Tin in alloys; Zinc; Zinc in insecticides; Hydrazine; Hydrogen peroxide; Peroxide in litharge; Other peroxides, dichromates, etc.; Removal of chloroform soluble colored substances; The determination of sodium thiosulphate; Tetrathionates; Sulphurous acids and sulphites; Bibliography; A partial list of atomic weights.

VOLUMETRIC IODATE METHODS

George S. Jamieson

96 Pages \$2.50

POLAR MOLECULES

P. Debye

172 Pages \$3.50

A survey of the notion of polar molecules as being systems having a distribution of electrical charges which can be characterized by a permanent electric moment. Absolute measurements of polarity may be used to determine the geometrical arrangements of the atoms in the molecules and by the quantum theory it is possible to develop the connection between polarity and the phenomena of dispersion and absorption especially in the infrared.

Contents: Fundamental electrostatic field relations; Polarizability and molecular structure; Measurements of polarity and its connection with chemical structure; Constitution of simple polar molecules; Anomalous dispersion for radio frequencies; Electrical saturation effects; Dielectric phenomena and quantum theory; Energy levels and wave mechanics; Rotating molecules; Dispersion and absorption of polar gases.

THE CHEMICAL CATALOG COMPANY, INC.

419 Fourth Avenue

New York, U. S. A.