

SCIENCE

7 December 1956

Volume 124, Number 3232

Editorial	Higher Education	1123
Articles	Scientific Research in the Soviet Union: <i>R. E. Marshak</i>	1125
	AAAS Meeting, New York: <i>R. L. Taylor</i>	1127
News of Science	Some Highlights of Astronomy in 1956; Fallout Hearings Proposed; Bioflavonoids and Colds; First AEC Awards Ceremony; New Neurological Research Foundation; Deep-Sea Diving Record; British Antarctic Expeditions; Scientists in the News; Recent Deaths; Education; Grants, Fellowships, and Awards; In the Laboratories; Miscellaneous	1140
Reports	Graying of Hair in Rats Fed a Ration Deficient in Lysine: <i>P. Vohra and F. H. Kratzer</i>	1145
	Spontaneous Rewarming of the Hypothermic Curarized Dog: <i>A. Y. Werner, D. Dawson, E. Hardenbergh</i>	1145
	Contamination of Human Cell Cultures by Pleuropneumonia-like Organisms: <i>L. B. Robinson, R. H. Wichelhausen, B. Roizman</i>	1147
	Comparative Studies of Lipoproteins by Starch and Paper Electrophoresis: <i>F. Paronetto, C. I. Wang, D. Adlersberg</i>	1148
	Missing Step in Guinea Pigs Required for the Biosynthesis of L-Ascorbic Acid: <i>J. J. Burns, P. Peyser, A. Moltz</i>	1148
	Uric Acid Metabolism and the Mechanism of Iron Release from Hepatic Ferritin: <i>S. Green and A. Mazur</i>	1149
	Criteria for Assessing Effects of Drugs on Posttetanic Potentiation: <i>D. W. Esplin</i>	1150
Book Reviews	<i>Atomic Quest; Kernmomente; Automation; New Concepts in Flowering-Plant Taxonomy; Modern Views on the Secretion of Urine</i> ; Books Reviewed in <i>The Scientific Monthly</i> , December; New Books; Miscellaneous Publications	1152
Letters	Quality and Quantity: <i>J. W. Still</i>	1155
	Research and Freedom: <i>W. S. Jarnagin</i>	1155
	Terminology: <i>D. L. Jofte and S. Warren</i>	1155
	Training of High-School Teachers: <i>E. Bozler</i>	1156
	Hidden History of a Dictionary: <i>W. C. Michels</i>	1156
	Australia Free from Rabies: <i>R. N. Wardle</i>	1156
Meetings and Societies	Anthropological and Ethnological Sciences; Social Sciences and International Peaceful Cooperation; AAAS 50-Year Members; Meeting Notes; Forthcoming Events	1157
	Equipment News	1169

Podbielniak CHROMACONS

Analyze Almost Anything Accurately


from HYDROGEN to LUBRICATING OILS

Fourteen Models


TO COVER EVERY SPECIFIC REQUIREMENT

SERIES 9485 (TWO MODELS)

Model 9485-5—Five Streams
Model 9485-1—One Stream
CHROMACON Series 9485, available in two models, designed for continuous (cyclic) analysis of gaseous or liquid process streams. Sample automatically charged, controlled by timer to permit any frequency of analysis from 1 to 12 analyses per hour. Operating temperature range 25° C. to 150° C., complete with recorder, detector cell and automatic charging device.

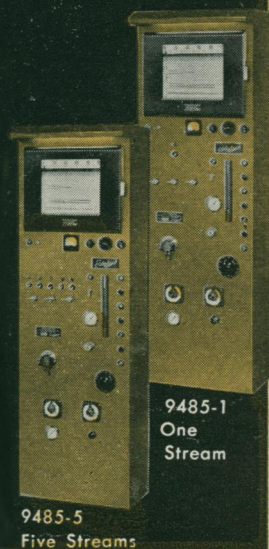

SERIES 9475 (SIX MODELS)

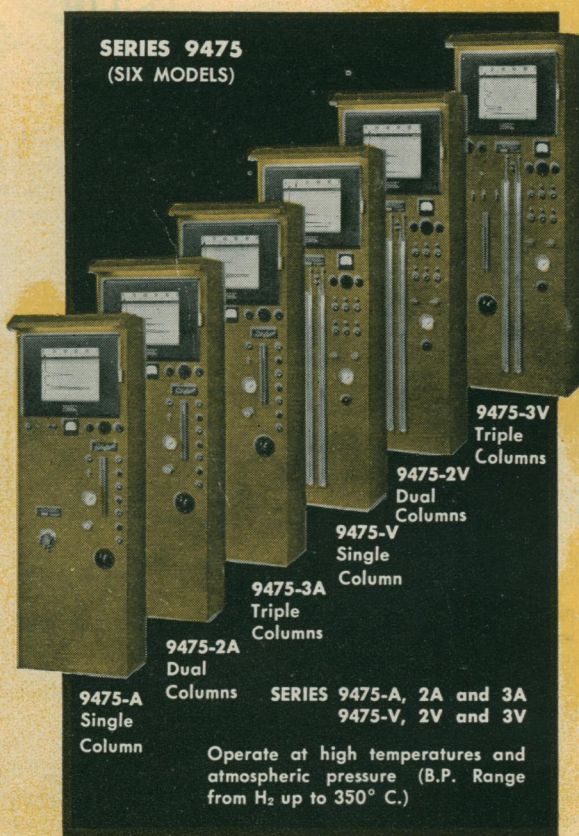

9475-A
Single
Column

9475-2A
Dual
Columns

9475-3A
Triple
Columns

9475-V
Single
Column


9475-2V
Dual
Columns

9475-3V
Triple
Columns

SERIES 9475-A, 2A and 3A
9475-V, 2V and 3V

Operate at high temperatures and atmospheric pressure (B.P. Range from H₂ up to 350° C.)

SERIES 9400 (SIX MODELS)

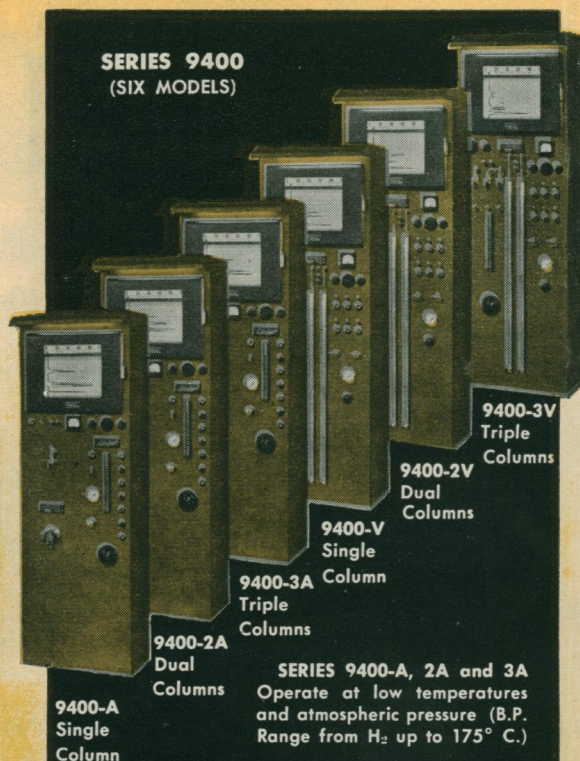

9400-A
Single
Column

9400-2A
Dual
Columns

9400-3A
Triple
Columns

9400-V
Single
Column

9400-2V
Dual
Columns

9400-3V
Triple
Columns

SERIES 9400-A, 2A and 3A
Operate at low temperatures and atmospheric pressure (B.P. Range from H₂ up to 175° C.)

SERIES 9400-V, 2V and 3V
Operate at low temperatures and subatmospheric pressure (B.P. Range from H₂ up to 200° C.)

COMMERCIAL PURIFICATION OF RARE PRODUCTS

Podbielniak Vapor-Phase CHROMANEER PLANT—Is used for the industrial purification of liquid or gaseous materials produced in small quantities. This method of separation gives a product of extreme purity, regardless of the boiling point of contaminating materials, with continuous automatic operation.

APPLICATIONS


1. Purification of silicone tetrachloride preparatory to the production of silicone transistors.
2. Production of flavoring oils, extracts, etc.
3. Production of perfumes.
4. Separation of pure compounds as in the production of "fine" chemicals or highly pure materials used as standards. Applicable to practically any organic liquid or gaseous material.

SPECIFICATIONS

1. COLUMNS: 2" to 4" diameter, 20 feet and up in length.
2. CAPACITY: 10 to 25 cubic feet of a gaseous product per day. Capacity varies with purity of charge and length of time cycle.
3. TEMPERATURE OF OPERATION: Columns thermostatted for operation at any temperature from 25° C to 200° C, as required for the particular separation to be made.
4. INSTRUMENTATION: Composition of exit gas monitored and continuously recorded.

WRITE FOR BULLETIN A-15 FOR COMPLETE SPECIFICATIONS.

CH-F-856 Copyright 1956 by Podbielniak, Inc. *Patents and Patents pending

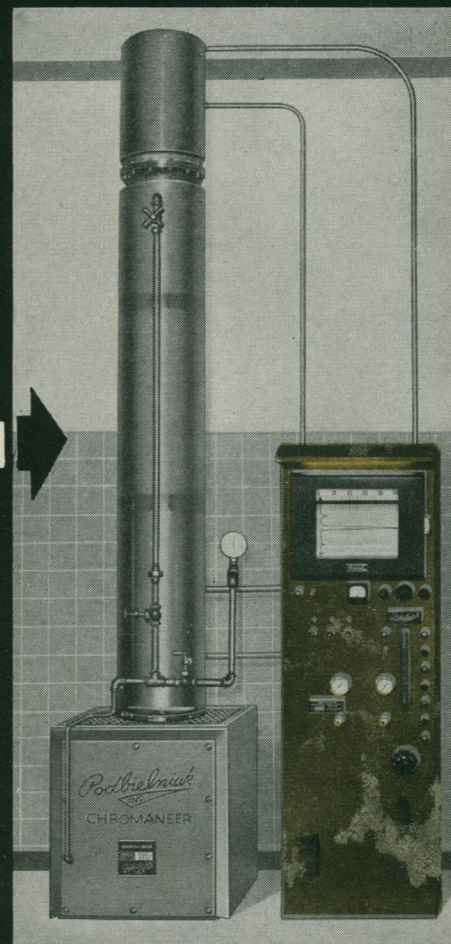


PODBIELNIAK, INC.


• 341 East Ohio Street •

Chicago 11, Illinois


Professional Illuminators

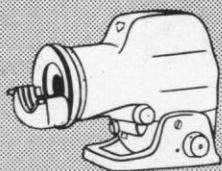


Model PR-27

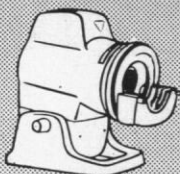

Model PG-26

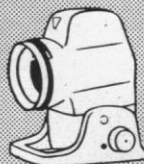

Model PL-25

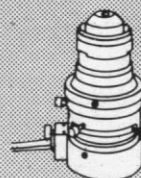

Panfocal
Illuminator


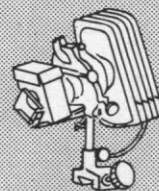



Photo-Micro
Illuminator



Grating
Monochromator

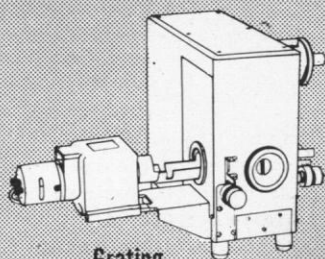


Nicholas
Illuminator

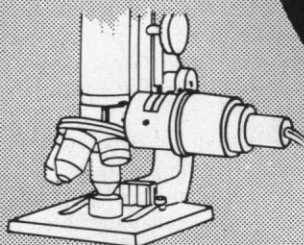


Tri-Vert Illuminator



Fluorescent
Illuminator

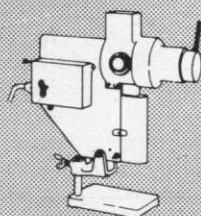


Arc Illuminator

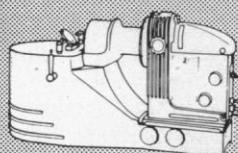


SpeedMatic Micro-projector
Basic Illuminating Unit

**TIME-SAVERS!
EFFORT-SAVERS!
SIGHT-SAVERS!**


**BAUSCH & LOMB
MICROSCOPE
ILLUMINATORS**

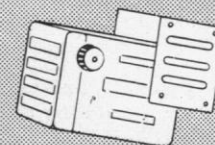


Opti-Lume
Illuminator

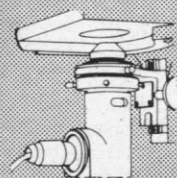

Dark Field
Illuminator

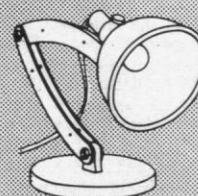

Reflector
Illuminator

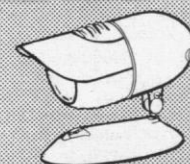

Micro-Lite

For vivid, detailed images . . . for comfort and convenience . . . choose the right light for your specific need. Choose from the most complete line serving science. Dependable illuminators for every visual and photographic application . . . for every budget.

WRITE FOR CATALOG DATA

Also, on request, obligation-free survey and evaluation of your present micro-illumination methods. Write to Bausch & Lomb Optical Co., 64212 St. Paul St. Rochester 2, New York.

BAUSCH & LOMB



America's only complete optical source . . . from glass to finished product.