

now you can weigh

MILLIGRAMS
in
SECONDS

to
MICROGRAM
ACCURACY
with the

CAHN ELECTROBALANCE



A major advance in weighing instrumentation combining a novel application of physical principles with modern engineering design, the Cahn Electrobalance is *unaffected by vibration, or temperature change*. You can use it anywhere, and depend on it.

The beam balances in a few seconds

Especially fast for weighing out pre-determined amounts.

One-Knob, one-dial operation

makes it exceptionally easy to use.

Precise and Accurate

It is really six balances in one.

| Range | Precision | Accuracy |
|-----------|---------------|------------------|
| 0-1 mg. | 1 microgram | 1.0 microgram |
| 0-5 mg. | 1 microgram | 2.5 micrograms |
| 0-10 mg. | 2 micrograms | 5.0 micrograms |
| 0-20 mg. | 5 micrograms | 12.0 micrograms |
| 0-50 mg. | 8 micrograms | 30.0 micrograms |
| 0-100 mg. | 25 micrograms | 100.0 micrograms |

Get complete information on this valuable new research tool by writing

CAHN Instrument Company

14511 Paramount Blvd., Paramount, Calif.

let us know if you would like a demonstration in your lab

new developments • new materials • integrated presentation of aliphatic and aromatic chemistry • excellent introduction to concepts of pi-bonding and conformation
in the 2nd edition of

A SHORT COURSE IN ORGANIC CHEMISTRY

Harold Hart and Robert D. Schuetz

approx. 344 pages

April, 1959

emphasis on principles and scientific thought • modern electronic approach to chemistry • modern concepts of energy science • influence of science on philosophic thought

in the 3rd edition of

THE STUDY OF THE PHYSICAL WORLD

Nicholas D. Cheronis, James B. Parsons, and Conrad E. Ronneberg

684 pages

1958

\$7.50

phylum-by-phylum approach • introduction to all biological concepts that have an important bearing on zoology • careful development of scientific ideas and experimental methods

in

GENERAL ZOOLOGY

Gairdner B. Moment

632 pages

1958

\$7.50

Laboratory Manual to accompany GENERAL ZOOLOGY prepared by H. Bentley Glass, Gairdner B. Moment, and Neal B. Weber

Available April, 1959

latest discoveries in the genetics of bacteria and viruses • new material on the nature of the gene • expanded section on the effects of radiation on genetics • new material on human heredity

in the 2nd edition of

GENETICS

A. M. Winchester

414 pages

1958

\$6.25

greater over-all clarity • logical organization of material • clear explanation of each step • excellent preparation for analytic geometry and calculus • new problems

in the 2nd edition of

PLANE TRIGONOMETRY

John J. Corliss and Winifred V. Berglund

397 pages

1958

\$4.00

HOUGHTON

Boston 7
Massachusetts

MIFFLIN

COMPANY



New York 16
New York
Atlanta 5
Georgia
Geneva
Illinois
Dallas 1
Texas
Palo Alto
California

Up IN THE
AIR
about
LABORATORY
COSTS?

Save:
TIME • SPACE • MONEY
WITH
ACE "MINI-LAB"

The Most Popular Line of Custom-Engineered, Interchangeable Small-Scale Glassware in existence.

Many items are equipped with $\$$ 18/9 Joints as well as the Popular $\$$ 14/20.

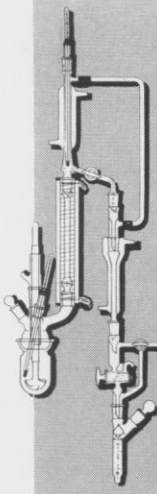Here is an Extremely Versatile Distillation Assembly consisting of several "Mini-Lab" Components. And—these same components lend themselves to a variety of uses. Ace stocks countless interchangeable assemblies and components. The individual user can readily develop innumerable variations of basic assemblies by using the many items available.

"Mini-Lab" uses only minimum amounts of your costly ingredients, and its compactness saves countless square feet of valuable lab space.

Ask Dept. ML-E to send you the
"Mini-Lab" Brochure.

ACE GLASS
INCORPORATED
VINELAND ♠ NEW JERSEY

Midwestern Division
LOUISVILLE, KY.—Box 996

New Products

The information reported here is obtained from manufacturers and from other sources considered to be reliable, and it reflects the claims of the manufacturer or other source. Neither Science nor the writer assumes responsibility for the accuracy of the information. A coupon for use in making inquiries concerning the items listed is included in the postcard insert.

■ **DIGITAL VOLTMETER** reads positive and negative d-c voltages from 100 mv to 999 v with automatic selection of range and polarity. Voltages are displayed as three digits with automatic decimal-point placement. Accuracy is ± 0.2 percent of reading ± 1 count. Input impedance to d-c is 11 megohm. Sampling rate is adjustable from 20 per minute to 5 per second and can be externally controlled. (Hewlett-Packard Co., Dept. 738)

■ **COINCIDENCE UNIT** features direct drive from photomultipliers, input limiting, and ultimate resolution of 3 msec with one anticoincidence and three coincidence channels. The coincidence criterion is set by a variable discriminator in a slow channel. All bias and supply potentials are electronically regulated. (E-H Research Laboratories, Dept. 745)

■ **IMAGE-CONVERTER CAMERA** is designed for use in photographing transient phenomena. The camera is electrostatically focused and deflected. Exposure time ranges from 0.05 to 1 μ sec with frame spacing from 0.5 to 50 μ sec. Number of frames is selectable from 1 to 5. Resolution at the center of the field ranges from 8 lines per millimeter for exposures from 0.05 to 0.15 μ sec to 17 lines per millimeter for steady exposure. Maximum image size at the anode is 25 mm for one to three frames and 15 mm for five frames. Spectral response is maximum at 4000 \pm 1000 Å. (Librascope, Inc., Dept. 741)

■ **PRESSURE SWITCH** consists of a contactor manometer and a relay and power supply. Control sensitivity of 0.005-in.-H₂O is achieved by use of manometer fluid of specific gravity 1. Relays operate on 110 v a-c and have double-pole, double-throw switching action. Contact capacity is 5 amp at 115 v a-c. (Meriam Instrument Co., Dept. 746)

■ **BLOOD-PRESSURE MONITOR** continuously indicates and records systolic blood pressure. Pressure in a digit cuff is maintained at the systolic blood pressure by means of a servomechanism. Pressure up to 300 mm-Hg may be determined. This pressure is continuously indicated on a manometer and can be fed to a pressure recorder if a permanent record is desired. The cuff does not reach completely around the finger, and thus a path for

L.L.U. TRANSLATIONS BULLETIN. Issued monthly by the Lending Library Unit of the Department of Scientific and Industrial Research. The Bulletin is a modification of the Translated Contents List of Russian Periodicals, which has ceased publication; the L.L.U. Translations Bulletin will list all translations which are produced by the D.S.I.R. Lending Library Unit. It will also contain information about translations available from other organizations. Articles on Russian Science and Technology will appear from time to time.

Annual Subscription **\$9.54**

VISUAL PROBLEMS OF COLOUR. An international Symposium held at the National Physical Laboratory. — *Even though the subject is a very complicated one, these proceedings—indispensable to the specialist—should also be very valuable to the non-specialist. They provide a full and well-balanced body of authoritative information at the present critical stage of our understanding of colour vision.* — *The New Scientist*, Feb. 19, 1959.

2 volumes **\$7.56**

THE MICROBIOLOGY OF FISH AND MEAT CURING BRINES. *Proceedings of the second international symposium on Food Microbiology, held at Cambridge, United Kingdom, in April 1957, edited by B. P. Eddy* **\$8.10**

BRITISH INFORMATION SERVICES

45 Rockefeller Plaza
New York, N.Y.

PSYCHOPHARMACOLOGY

AAAS Symposium Volume

edited by
Nathan S. Kline

6" x 9", clothbound, 175 pp.,
bibliographies, index, 1956
**\$3.50 (\$3.00 for
cash orders by AAAS members)**

This volume consists of material presented at the first major conference on the remarkably successful use of new drugs such as chlorpromazine in the treatment of mental disease.

"The authors are competent and careful workers who have approached this problem with a scientific attitude . . . Throughout the volume runs the thread of caution . . . New vistas are being opened for the psychiatrist, the neurologist, the physiologist, the psychologist, the pharmacologist, and the chemist." Foreword, Winfred Overholser.

"This volume is not a reference intended for use at the introductory student level. It can be reviewed with interest, however, by any serious member of the reading public." *American Journal of Pharmaceutical Education*, July 1956.

At all book stores or write

**American Association
for the
Advancement of Science**
1515 Massachusetts Ave., NW
Washington 5, D.C.

return of blood to the heart is provided. Timers permit automatic control of duty cycle. The associated recorder is a pen-writing type with 5-in. chart width. Chart speeds are 0.2 and 1 in./min. (Medical Electronics Development Co., Dept. 749)

■ **ULTRASONIC IMPACT GRINDER** uses a magnetostrictive transducer to permit a 100-percent duty cycle. Materials that can be worked are semiconductors, ferrites, carbides, metals, jewels, and other hard and brittle materials. Cutting force is adjustable. The slurry tray is readily removable, and tool cones are easily interchanged. The ultrasonic driver delivers 100 w. (Raytheon Manufacturing Co., Dept. 750)


■ **SCALER SYSTEM** offers preset count from 20 to 10^7 and preset time from 0.02 to 1000 min with automatic operation. Scaling operations stop when either preset time or count is reached, whichever occurs first. The system automatically commences a new count after print-out. Any number of units can be used with a single printer; data are read out sequentially. (Radiation Instrument Development Laboratory, Inc., Dept. 755)

■ **CARBON ANALYZER** is semiautomatic, using a single source of oxygen to burn the sample and to raise and lower the concentration in solutions. The apparatus features replaceable float valves, screw-type zero adjustment, and built-in circulating pump. Operation is on 115-v, 60-cy/sec power. (Laboratory Equipment Corp., Dept. 748)

■ **COMPRESSED-GAS REGULATOR** features small size: 1 by $2\frac{3}{4}$ by $3\frac{1}{2}$ in. The regulator, a two-stage device, will handle oxygen, nitrogen, or other gases. Delivered-pressure range is 0 to 90 lb/in². (Beacon Devices, Inc., Dept. 747)

■ **VACUUM-GAGE CONTROL** provides readings from one ionization and two thermocouple gages covering the range from 1 to 10^{-8} mm-Hg. A protective relay guards the ionization gage against burn-out due to pressure surges and permits actuation of warning devices. Ionization emission currents of 1, 5, or 10 ma may be selected. (NRC Equipment Corp., Dept. 754)

■ **CURRENT INTEGRATOR** accepts input currents in six decade ranges from 0.001 to 100 μ a full scale. Integration range varies, with current range, from 1 μ coul to 10 coul with 100-coul range available on request. Accuracy of current measurement and of charge measurement is ± 0.05 percent of reading. Current is amplified by a feedback-stabilized current amplifier. A capacitor integrates the amplifier output to a preset voltage

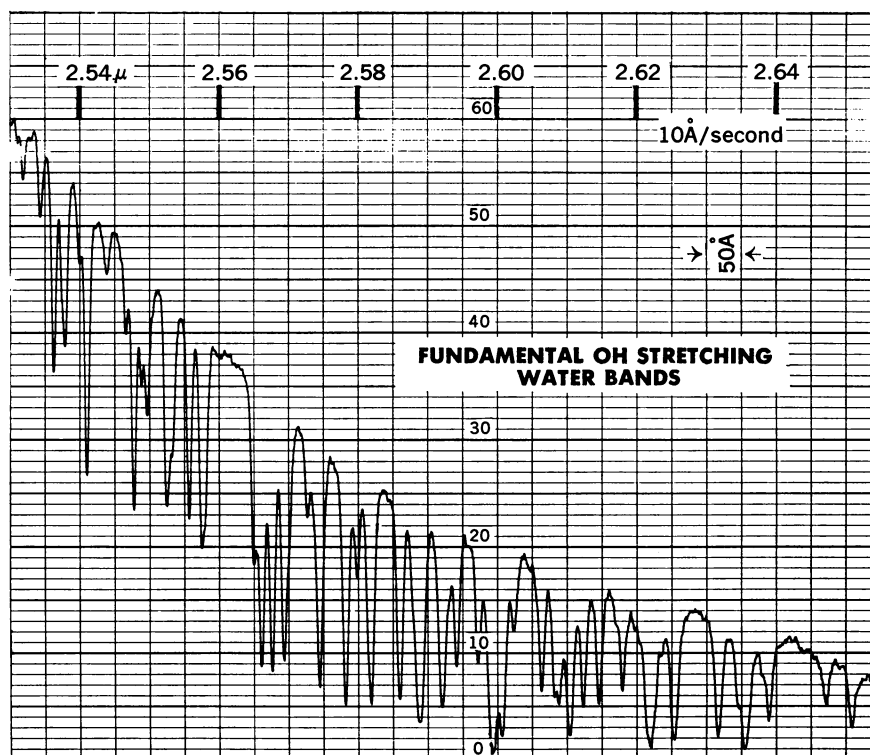


FUNDAMENTAL OH STRETCHING WATER BANDS

RANGE: 2.50-2.65 microns RESOLUTION: 1 cm^{-1} (6\AA)

The above curve illustrates the high resolution that users of Cary Model 14 Spectrophotometers are getting for measurements in the near-IR region. The Model 14's ability to resolve such fine structure is a feature not ordinarily found in a general purpose instrument having a wide wavelength range (1860 \AA -2.65 microns). In most of the ultraviolet-visible region, resolving power of the Model 14 is better than 1 \AA .

High resolving power is just one of many features that make the Model 14 so useful. A broad wavelength range, a wide choice of scanning and chart speeds, accommodation of a variety of types and sizes of sample cells, stray light of less than 1 ppm, photometric reproducibility better than .004 in absorbance even at high absorbance, and many special accessories suit the Model 14 to a wide variety of spectrophotometric problems requiring fast, accurate analyses. These and other performance features have made it the preferred recording spectrophotometer of leading research laboratories throughout the world.

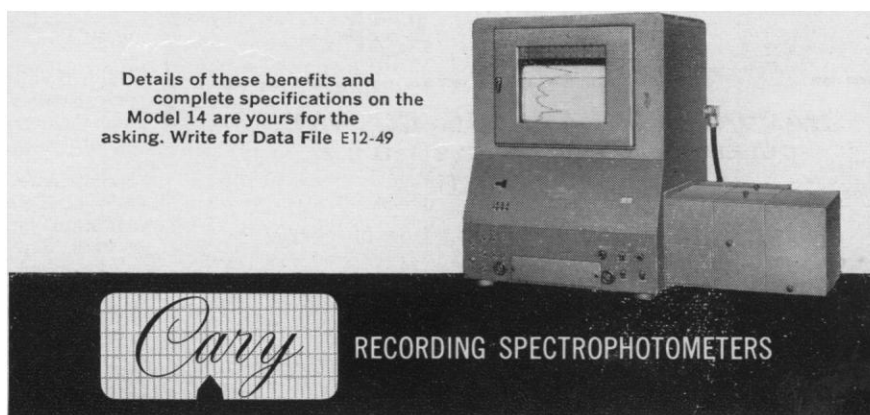

APPLIED PHYSICS CORPORATION • 2724 SOUTH PECK ROAD • MONROVIA, CALIFORNIA

at which it automatically discharges. The number of discharges and the final charge on the capacitor are recorded. An internal calibration signal, which is applied to the input of the current amplifiers tests for proper operation. (Applied Radiation Corp., Dept. 765)

■ **TRIPLE-POINT CELL** maintains a reference temperature of 0.01°C by establishing equilibrium between ice, water, and water vapor by partially freezing pure water in an evacuated-and-sealed glass system. Cells are furnished with full instructions for use. (Transonics Inc., Dept. 751)

■ **GAUSSMETER** is transistorized and designed to operate either on batteries or external power supply. Fields up to 30,000 gauss are measured with indication of direction. No ballistic reading is given. Probe is 0.025 in. thick and has an active area of 0.01 in². (Dyna-Empire, Inc., Dept. 762)

■ **VACUUM PUMP** for vacuum-distillation applications maintains a pressure of 0.015 mm-Hg with a free-air capacity of 36 lit./min at 540 rev/min. The open-chamber design eliminates freezing and permits solid matter to be exhausted. (Nelson Vacuum Pump Co., Dept. 760)

■ **CONTINUOUS-WRITING FRAMING CAMERA** operates with no blind time and thus can be used to photograph high-speed events that do not lend themselves to synchronization with the position of the rotating mirror in the camera. Magnification-velocity products up to 3 mm/ μsec are accommodated. At conventional magnifications, object velocity can be Mach 35. Daylight loading with standard 35-mm film cassettes is provided. (Beckman and Whitley, Inc., Dept. 758)

■ **WHITENESS REFLECTOMETER** uses paired vacuum tubes in a null bridge circuit. It employs 45 deg-0 deg geometry and has green and blue tristimulus filters, each with a separate pair of phototubes. Separable sections contain the measuring and indicating circuit and the sample viewing system. (Hunter Associates Laboratory Inc., Dept. 752)

■ **VIBRATION GENERATORS** are available in two models rated at 500 and 150 lb force, respectively. The smaller model allows a bare-table resonance of more than 7000 cy/sec. Bare-table resonance of the larger model is 6000 cy/sec. The generators are of electrodynamic design with dual magnetic field structure producing low stray fields. (Ling Electronics, Inc., Dept. 759)


■ **MICROWAVE-FREQUENCY CALIBRATOR** generates harmonics up to 25 kMcy/sec. The instrument's 450 Mcy/sec crystal-controlled signal is designed to be fed directly into a waveguide or coaxial crystal holder. A 5-Mcy/sec crystal permits calibration against station WWV. (Micro-Now Instrument Co., Dept. 763)

■ **INVERTERS**, using static elements, operate from 28-v d-c input and provide output of 115 v at 400 cy and 2000 cy/sec. Regulation characteristics are ± 1 percent for voltage and ± 0.02 percent for frequency. Distortion is 5 to 10 percent. Weight ranges from 5 to 15 lb. (Gulton Industries, Inc., Dept. 764)

■ **MILLIOHMETER** features 13 overlapping ranges measuring 0.001 to 1000 ohm full-scale. Accuracy is ± 3 percent of full scale except in the 1-milliohm range, where it is ± 5 percent. No calibration or zero adjustments are required. The instrument applies only 2 μw across the circuit under test. (Keithley Instruments Inc., Dept. 756)

■ **DIRECT-CURRENT METER** has the equivalent of a 63-in. scale with 1000 divisions. This scale is an expansion of the actual 6.3-in. scale effected by a circuit that permits the instrument to cover any 10

New CATALOG NOW AVAILABLE

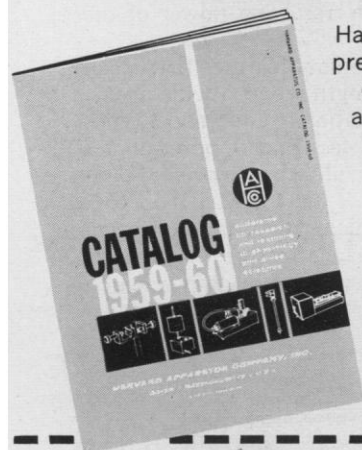

Harvard Apparatus Co., Inc. presents its 1959-60 catalog of apparatus for research and teaching in physiology and allied sciences. It features many new products as well as our complete line of equipment and supplies. The catalog is available free of charge.




HARVARD APPARATUS CO., INC.
DOVER • MASSACHUSETTS • U. S. A.
(a non-profit organization)

Gentlemen:
Please send new 1959-60 catalog free of charge.

Name _____
Title _____
Organization _____
Address _____
City _____ Zone _____ State _____

FERMENTOR DRIVE ASSEMBLY FOR LABORATORY AND PILOT PLANT APPLICATION

- Aerobic or Anaerobic Fermentations • Tissue Cultures
- Other Metabolic Studies in Submerged Cultures

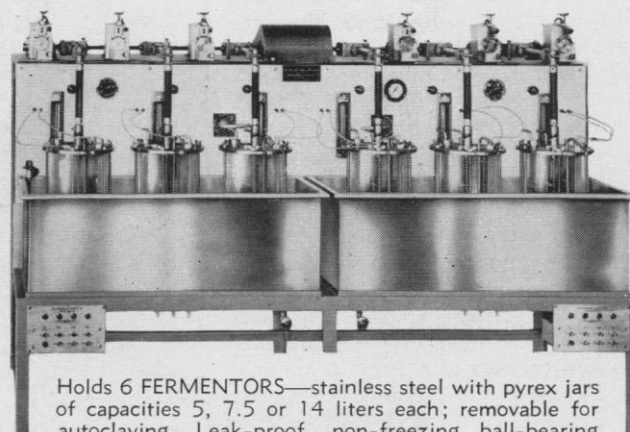

Holds 6 FERMENTORS—stainless steel with pyrex jars of capacities 5, 7.5 or 14 liters each; removable for autoclaving. Leak-proof, non-freezing ball-bearing agitator housings.

Stainless steel baths, thermostatically controlled within $\pm \frac{1}{2}^{\circ}\text{C}$. Twin, anti-friction drives for wide range of agitation rates. Available with automatic anti-foam system. Single fermentor and 3-fermentor assemblies also available.



NEW BRUNSWICK SCIENTIFIC CO.
PRECISION LABORATORY APPARATUS
P. O. BOX 606 • NEW BRUNSWICK, N. J.

Write for Bulletin F-S424



CHICAGO, U.S.A.

STUDENT MICROSCOPES

Most reasonably priced **GUARANTEED**
Microscope on the market.

Made in West Germany

NEW DESIGN
EXCLUSIVE
SAFETY FEATURES
HIGH QUALITY OPTICS
10X OCULAR
OBJECTIVES

16mm (10X) N.A. 0.27

4mm (44X) N.A. 0.66

plano/concave mirror

STILL \$118.00

TEN YEAR GUARANTEE

Write for catalogue
listing safety features


10% Discount on 5 or more.
Models may be assorted to ob-
tain this discount

TRANSPORTATION
INCLUDED

THE GRAF-APSCO CO.

5868 BROADWAY

CHICAGO 40, ILL.

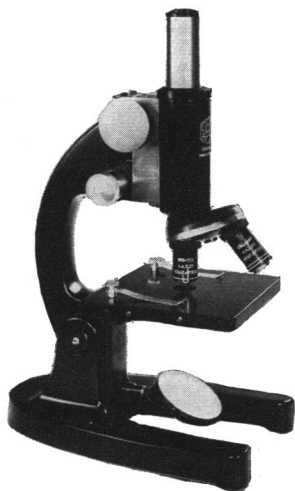

Accuracy to 1/10 Microliter!

HAMILTON MICROLITER SYRINGES WITH CEMENTED NEEDLES

When *accurate* liquid discharge is vital, such as in gas or vapor chromatography, these Hamilton Syringes are essential to your analytical work . . .

- * Three models—100, 50 and 10 microliter capacities—meet all requirements
- * Skilled operators maintain accuracy of discharge to 1/10 of 1 microliter
- * Precision of graduation is consistently better than 1%
- * Fixed concentric needle minimizes dead volume in the syringe, prevents leakage
- * Syringe is tested to be leak-tight with water at 150 psi
- * Special needle point design easily perforates rubber closures without plugging

Order direct—or write *today* for literature and prices!


THE HAMILTON COMPANY

DEPARTMENT K

1134 WHITLEY STREET • WHITTIER, CALIFORNIA

*A Complete Line of Precision,
Radiochemical Handling Equipment*


New COLOR-EYE MODEL D

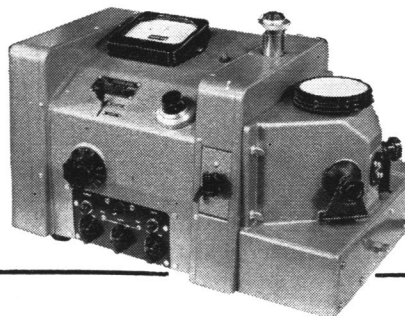

PRECISION COLORIMETER and ABRIDGED SPECTROPHOTOMETER

DELTA E COMPUTER


is a COLOR-EYE accessory which computes deviations from a color standard and provides a single-value color tolerance in NBS units.

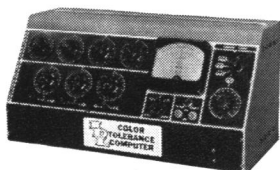

Write for technical literature.

- ▶ COLOR-EYE is a dual-purpose precision instrument for establishing and maintaining quality color control under laboratory and production conditions.
- ▶ The TRI-STIMULUS SYSTEM measures color differences on the three-color basis established by the C.I.E. It quickly determines the color differences in hue, value and chroma. Enables accurate production inspection and precise checking on color standards.
- ▶ The ABRIDGED SPECTROPHOTOMETER employs ten filters for analyzing color formulations and determining metameric conditions. It permits matching of colors, even with different types of pigment and different materials . . . textiles, plastics, painted surfaces, etc.
- ▶ A built-in COOLING SYSTEM controls sample temperature rise so that rapid analysis of colors in temperature-sensitive materials is possible.
- ▶ A replacement-type FILTER cleans the air to prevent dust from collecting on the optics and enables accurate measurement in dusty atmospheres.

*Demonstrations of Color-Eye in operation in
your plant can be arranged upon request.*

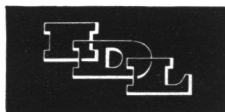

INSTRUMENT DEVELOPMENT LABORATORIES, INC.

An Affiliate of the Royal McBee Corporation

67 Mechanic Street, Attleboro, Massachusetts, U. S. A.

NEW liquefied gas container KEEPS CONTENTS TWICE AS LONG!



Only 22 $\frac{3}{4}$ in. high, 15 $\frac{3}{8}$ in. diam., tough, corrosion-resistant welded aluminum construction cuts weight to 19 lbs. Low conductivity neck tube has 1 $\frac{1}{2}$ in. opening.


Compared to ordinary vacuum-insulated vessels, the LINDE LD-25 has a 50% lower evaporation rate! It's the most practical and economical container available, and can bring substantial savings in the storage of liquefied atmospheric gases. Here's why:

- Holding time—liquid oxygen loss rate at only 4% per day; nitrogen at 5%; neon, 4%; argon, 8%.
- Capacity is 25 liters, yet weighs only 19 lbs. empty—less than containers that hold only $\frac{1}{2}$ as much.
- Wide neck permits emptying in 90 seconds, 10 times faster than old-fashioned narrow-neck containers. Filling time is also reduced.
- Automatic pressure withdrawal tube, dipper, and roller caster base are available.



LINDE guarantees all of its containers against defective material and workmanship for a period of one year from date of shipment. This includes a one year guarantee against excessive evaporation loss.

Linde

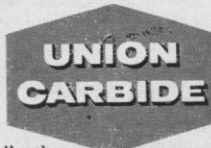

"Linde" and "Union Carbide" are registered trade-marks of Union Carbide Corporation.

Linde Company, Division of Union Carbide Corporation Dept. SC-44
30 East 42nd Street, New York 17, N.Y.

Please send me information on

- ☐ LD-25 Liquefied Gas Container
☐ other equipment for liquefied atmospheric gases
☐ (please specify) _____

Name _____

Firm Name _____

Address _____

City _____ Zone _____ State _____

percent of full-scale range over a full 6.3-in. scale. The meter is direct-reading. Readability and repeatability are ± 0.01 percent. Accuracy is ± 0.05 percent of full-scale range. Single-range or switch-controlled multirange models are available. (Sensitive Research Instrument Corp., Dept. 757)

■ HALL GENERATORS allowing multiplication of two electrical quantities feature better than 1-percent linearity over an output range of 100,000 to 1. Units, made of indium arsenide and indium antimonide, are available with external magnetic field coils or in the form of probes. (Device Development Corporation, Dept. 753)

■ TEMPERATURE RECORDER provides a variety of ranges by means of plug-in elements matched to appropriate thermocouples. Automatic reference junction compensation eliminates the need for external temperature reference. Accuracy is ± 1 percent of range. Chart width is 5-in., and 1-, 2- and 4-chart-speed models are available. (Varian Associates, Dept. 761)

■ SHAFT-ROTATION-TO-DIGITAL CONVERTER is available in ranges 0 to 999 or 0 to 1999. Coding is accomplished by a disk with read-out brushes not normally in contact. Read-out command places brushes in contact and then energizes the disk. Ambiguity is prevented by a special mechanism. Any type of code can be provided. (Norwood Controls, Dept. 768)

■ HANDLING TOOL for delicate parts is vacuum-operated. A solenoid-operated valve between the vacuum line and the pickup is controlled by a foot-operated switch, by which the operator can release the suction when the part is correctly positioned. Operation is on 115-v, 60-cy/sec power. (Carmen Laboratories, Dept. 766)

■ PULSE GENERATOR produces pulses of controlled amplitude and duration at frequencies from 10 cy to 5.0 Mcy/sec corresponding to the frequency of the input signal. Pulse duration is variable from 0.1 to 5 μ sec; rise time is 0.04 μ sec and fall time is 0.05 μ sec at 40 v output. Amplitude is 2 to 4 v across a 1000-ohm load. Input amplitude must be 3 v r.m.s. or greater. (Technitrol Engineering Company, Dept. 776)

■ GAS CELL for infrared absorption spectrophotometry employs a folded optical path 100 cm long. Sample volume is 950 cm³. The cell is constructed of cast aluminum with Neoprene gaskets and can contain pressures above or below atmospheric. (Perkin Elmer Corp., Dept. 779)

■ **TEMPERATURE CONTROLLER** with thermocouple input has a set-point resolution of 0.1 percent of range and handles up to 330 kva through standard saturable reactors. A range selector allows manual switching for temperatures 0° to 800°, 800° to 1800°, and 1800° to 2800°F. In a typical application the controller is said to limit fluctuations to 0.25°F at a temperature of 2200°F. (Hagan Chemicals and Controls Inc., Dept. 777)

■ **OPTICAL VIEWER** for observation of inaccessible regions features a 180-deg hemispheric field. Lengths up to 40 ft are available. Underwater design, sectional construction, special illumination, and photographic attachments can be provided. Depth of field extending from ¼ in. to infinity is claimed. (Lerma Engineering Corp., Dept. 769)

■ **SPECTRUM ANALYZER** features r.m.s., average, and peak read-outs. Center frequencies of the analyzer's 30 one-third octave filters range from 40 to 32,000 cy/sec. Filter characteristics are: tops flat within ±0.5 db; sides with maximum slope 120 db per octave; skirt selectivity greater than 50 db per octave and greater than 70 db two octaves from the center frequency. (B & K Instruments, Inc., Dept. 774)

■ **PIPETTE PUMP** accommodates all sizes of pipettes in rubber-lined chuck jaws. In addition to the large rubber bulb used for filling, a small internal bulb, depressed by a control-knob screw, provides control of micro quantities. The device can be dismantled for cleaning and sterilizing. (Arthur S. LaPine and Co., Dept. 767)

■ **SHOCK TESTER** develops thrust through differential gas pressure acting on two faces of a thrust piston in a closed cylinder. The cylinder is separated into chambers by an orifice plate. To control the thrust developed, a metering pin attached to the underside of the thrust piston projects through the orifice and meters the flow of gas as a function of piston motion. Maximum thrust of a 3-in. tester is 10,000 lb. (Consolidated Electrodynamics, Dept. 778)

■ **STREAKING PIPETTE** for setting up paper chromatograms is designed to eliminate drop formation at the tip of the pipette and to maintain a constant hydrostatic head during delivery of the entire sample. Approximately 0.3 ml may be applied by the pipette with one filling. (Research Specialties Co., Dept. 782)

■ **RESISTANCE SET** consists of 21 resistors scaled in a modified binary-decimal sequence of 1,2,2,4 by decades with an extra units' digit, permitting synthesis of

RECENT PERGAMON PUBLICATIONS

CERENKOV RADIATION and its Applications

by J. V. JELLEY
An up to date review of the theory, experiment and application of the phenomena connected with this form of radiation. With bibliography. Published for the United Kingdom Atomic Energy Authority. 304 pages illustrated \$10.00

CARBON: Proceedings of the Third Conference

Edited by S. MROZOWSKI
This volume covers the conference held at Buffalo in 1957, sponsored by the University of Buffalo, the National Science Foundation and the U.S. Office of Naval Research. 300 pages \$20.00

BORON FLUORIDE

And Its Compounds as Catalysts in Organic Chemistry

by A. V. TOPCHIEV, S. V. ZAVGORODNII, Y. M. PAUSKHIN
members of the USSR Academy of Sciences, Petroleum Institute.
Translated by J. T. GREAVES

The most complete and up to date study of the use of boron fluoride compounds as catalysts and of their derivatives and molecular compounds. 450 pages illustrated \$12.00

VISTAS IN FREE RADICAL CHEMISTRY

A Volume in Memory of the late Morris S. Kharasch. Edited by W. A. WATERS
This volume has permanent value as an accessible historical record of one of the major achievements of twentieth century organic chemistry. (Supplement Number 3 of 'Tetrahedron') 200 pages illustrated \$12.50

LIQUID SCINTILLATION COUNTING

Edited by C. G. BELL
A conference held at Northwestern University in 1957. It is the first volume on liquid scintillation counting, its theory, techniques and applications. 304 pages \$10.00

THEORY OF RELATIVITY

by W. PAULI, translated by G. FIELD, University of Birmingham
This book gave a complete review of the whole literature on relativity theory when it was originally published in 1921. It has been brought up to date by the late Professor Pauli for this English edition by adding 27 Supplementary Notes, giving selected information about later developments and his own view on some controversial questions. 242 pages \$6.00

PROGRESS IN METAL PHYSICS VII

Edited by B. CHALMERS
Volume 7 is now ready. The annual appearance of a new volume, with authoritative reviews of the existing state of knowledge, is keenly anticipated by a wide circle of faithful readers. 400 pages \$16.00

SCIENTIFIC MANPOWER IN EUROPE

by E. MCCRENSKY
This most clearly written and important book, is a comparative study of Scientific Manpower in the public service of Great Britain and selected European countries. It is the first book to describe the methods to recruit, retain, train and manage technicians and scientists for the public service. 188 pages \$6.50

SOVIET REVIEWS OF NUCLEAR SCIENCE


Being a translation of articles published in the Soviet journal *Atomnaya Energiya* commemorating the 40th Anniversary of the October Revolution. It is written by leading scientists who discuss the subject from the practical and theoretical angles. (International Series of Monographs on Nuclear Energy, Div. 11, Vol. 3) 110 pages \$5.00

PERGAMON PRESS INCORPORATED

122 East 55th Street, New York 22, N.Y.

WHAT'S BACK OF THE DK-2* ?

*BECKMAN
RATIO RECORDING
SPECTROPHOTOMETER

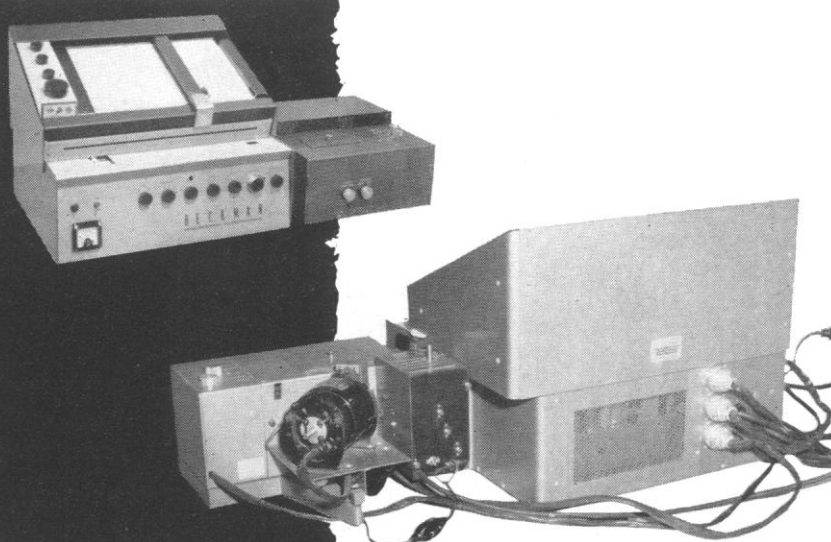

You, the alert scientist, are aware of the many specific features . . . all the time saving advantages . . . and the new capacity for your laboratory that are offered by the Beckman DK-2 Spectrophotometer.

Yes, the DK-2 is a new and complete 3 in 1 spectrochemical lab in itself . . . embracing UV . . . VISIBLE . . . and the near INFRARED.

But, do you know all that's back of the DK-2 installed by Will?

The Plus Values of a Will-Installed DK-2

With every Will-supplied DK-2 you receive these added values and services at no extra cost:

- 1st PLUS—Installation by a factory trained specialist, including performance checks and wavelength calibration.
- 2nd PLUS—Service for one year by Will's factory trained, Beckman-certified technicians.
- 3rd PLUS—Free inspection of your DK should you suspect defective parts or inferior performance.
- 4th PLUS—Prompt replacement of faulty materials.
- 5th PLUS—A free six months' subscription to the Anderson Laboratory Near-Infrared Spectra Program (value, \$150).

IMMEDIATELY AVAILABLE FROM WILL'S STOCK

- 24128—Beckman DK-2 Ratio Recording Spectrophotometer, complete with flat bed recorder, power supplies, dual source unit, detectors, matched silica cells and assorted chart paper; for 115 Volts, 60 cycle AC \$7975
- 24134—Beckman DK-2, same as above with addition of Linear Wavelength Presentation \$8850

To be sure you receive all that is back of the DK-2 . . . order from Will!

A progressive past . . .



A progressive future

Will CORPORATION
and subsidiaries

Specialists in  Scientific Supply

ROCHESTER 3, N.Y. • ATLANTA 1, GA. • NEW YORK 52, N.Y. • BALTIMORE 24, MD.
BUFFALO 5, N.Y. • SO. CHARLESTON 3, W. VA.

resistances from 0 to 100,000 ohm in 1-ohm steps. Accuracy and stability of ± 0.0015 percent of the maximum value are claimed. An oil bath is used to minimize thermal and self-heating effects. The assembly is housed in a hermetically sealed metal case. (Julie Research Laboratories, Inc., Dept. 780)

■ **FLUOROMETER** offers sensitivity to as little as 5×10^{-4} μg of quinine sulfate. A range multiplier covers 1:3:10:30 concentration ranges. Neutral-density filters provide additional range multiplication. A single multiplier phototube is used with mechanically chopped radiation of either 2537 or 3600 Å. (G. K. Turner Associates, Dept. 770)

■ **D-C AMPLIFIER**, fully transistorized, has gain of 10 to 10,000. Band width is from d-c to 60,000 cy/sec. Noise level is 6 μv . The amplifier may be used in single-ended or balanced configurations. (Princeton Science Associates, Inc., Dept. 771)

■ **THERMOELECTRIC GENERATOR** operates with efficiency of about 6 percent. Semiconducting materials are used in the thermoelectric elements. A 5-watt size is intended for use with radioisotope heat source and air cooling. (Minnesota Mining and Manufacturing Co., Dept. 772)

■ **CHARTING TOOL** is used to place narrow tapes (1/64 to 1/8 in. wide) on charts or graphs. A cutting bar permits the tape to be cut at any point. Tapes are available in 14 solid colors and ten distinct patterns. (Labelon Tape Co., Dept. 773)

■ **DIGITAL MILLIVOLTMETER** is a self-balancing potentiometer instrument with accuracy of ± 0.1 percent. Sensitivity of the basic model is 2.5 μv . Full-scale ranges from 10 mv to 100 v are available. Output of the potentiometer is converted to a three-decade digital indication for continuous visual display. The instrument can be equipped to operate printers, punches, or electric typewriters. Print-out can be made on command, or an optional timer can be supplied to actuate periodic print-out. (Computer Equipment Co., Dept. 781)

■ **VACUUM FURNACE** reaches temperatures up to 2700°F. Cylindrical radiation shields of heat-resisting alloy are arranged in several concentric shells to insulate the inside walls. Top and bottom have horizontal alloy sheets to form a complete enclosure around material being treated. The furnace casing is water jacketed. (General Electric Co., Dept. 798)

JOSHUA STERN
National Bureau of Standards,
Washington, D.C.

*announcing
the*

SERVALL®

*The latest word
in automatically-
controlled,
low-temperature
centrifuges*

RC-2



- **AUTOMATIC CONTROLS**
- **SUPERSPEEDS**
- **RCF IN EXCESS OF 35,000 x G**
- **CONTINUOUS FLOW**
- **MODERN, FUNCTIONAL DESIGN**

No fewer than five different rotors may be used in the new SERVALL RC-2: Superspeed, Large Capacity, Multi-Tube Superspeed, Horizontal, Virus & Particle Counting (with others under development). Gyro-Action Direct Drive: the only significant self-centering development in a decade; provides smoother operation than any other drive system. Continuous Flow System permits gallon quantities of material to be separated *directly* in tubes. Automatic Acceleration, Dynamic Braking, and exclusive Dual Automatic Temperature Control for accurately maintaining material temperature within $\pm 1^\circ\text{C}$ at all times.

Write for more information about this, the Superspeed Automatic Refrigerated Centrifuge that researchers everywhere have been asking for: Bulletin SC-4RC-2.



An independent company;
not connected with any
other centrifuge manufac-
turer. Established 1934.

Ivan Sorvall, Inc.


NORWALK • CONNECTICUT

DESIGNERS, MANUFACTURERS AND
DISTRIBUTORS OF "SERVALL"
LABORATORY INSTRUMENTS


Here's why AO Phasestar offers you more!

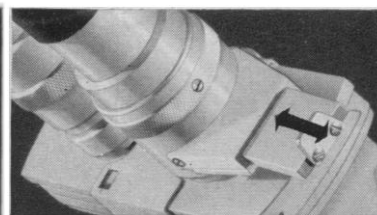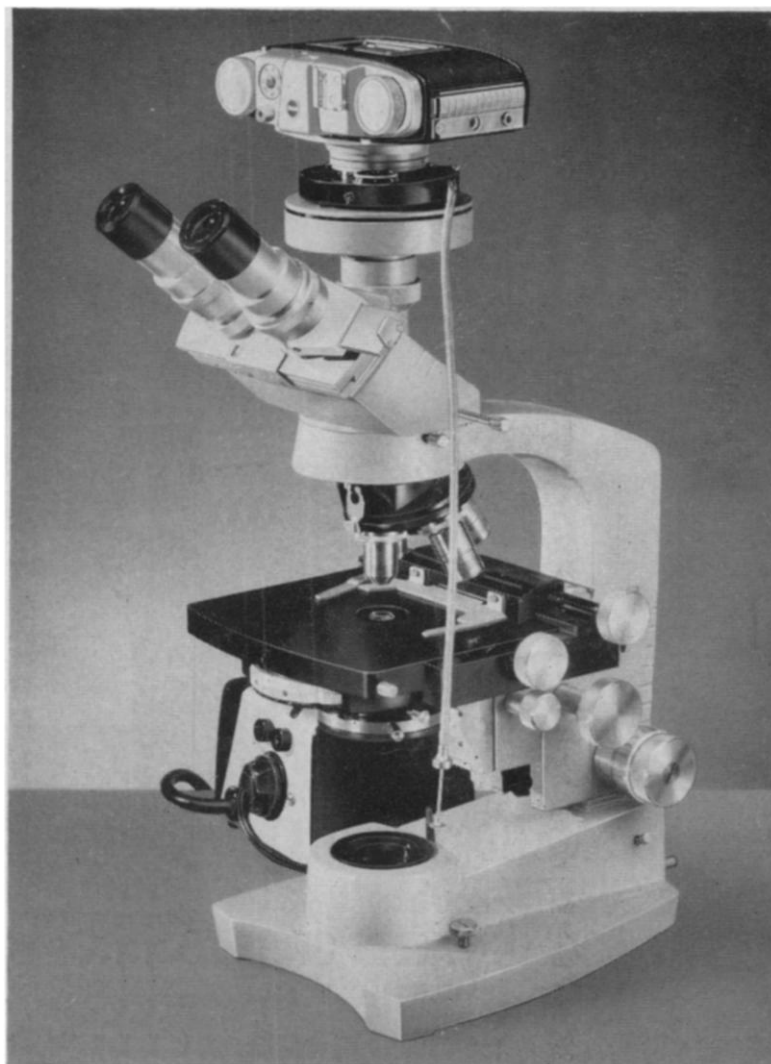

New Bertrand Lens Built-in Bertrand lens is now supplied as standard on all binocular and trinocular Phasestars to provide new convenience. You slide lens into position in eyepiece tube and raise or lower to focus quickly on diffraction plate at back aperture of objective.

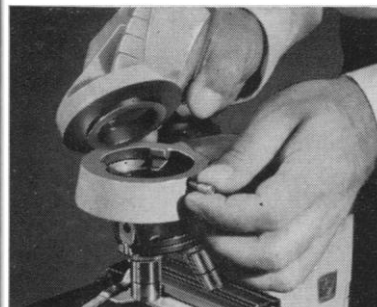

Interchangeable, Rotatable Bodies Your choice of interchangeable, rotatable, monocular, binocular or trinocular bodies. Simply loosen thumb screw to rotate (full 360°) or interchange. Special inclined monocular body with vertical photographic tube also available.

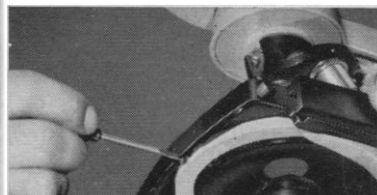

Interchangeable Stages You loosen single set screw to remove and interchange stages... your choice of graduated or ungraduated mechanical or MICRO-GLIDE stages.


The AO Spencer Phasestar is the most versatile phase microscope ever developed. Beginning with a basic Phasestar stand, you can select from a wide choice of interchangeable bodies, stages, optics and phase accessories to make up the specific model that meets your exact needs...all in all, more than 50 possible combinations or models to choose from. And whatever the model you choose, you're assured of superior AO Spencer quality *plus* matchless convenience, comfort and versatility.

Often particular phase applications require specific equipment and phase contrast objectives. Before you select your phase microscope and accessories consult your AO Sales Representative. He's qualified to help

you select the microscope and equipment you need to do the job. There's no obligation, of course.

AO SPENCER PHASE ACCESSORIES FOR PLATELET COUNTS

The use of phase contrast for blood platelet counts is increasing in acceptance and application. The longer working distance required with the phase hemacytometer can be obtained with special AO Spencer long working distance accessories that provide up to 7mm working distance in air or 10½mm in glass.

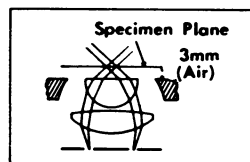

**American Optical
Company**

INSTRUMENT DIVISION, BUFFALO 15, NEW YORK

Dept. P-4

- ☐ I would like to review a specific phase microscopy problem. Please have AO Sales Representative call.
- ☐ Please send copy of Phasestar Brochure.

Name _____

Address _____

City _____ Zone _____ State _____