

28 May 1982 • Vol. 216 • No. 4549

\$2.50

SCIENCE

AMERICAN ASSOCIATION FOR THE ADVANCEMENT OF SCIENCE



The Gentle, Rapid Way to Separate Living Cells.

This technician has just prepared a fraction containing 95% granulocytes from human serum with the Beckman Elutriator Rotor System. Cell viability was excellent.

In the past four years, the technique of centrifugal elutriation has become extremely popular for synchronization studies of cultured cells, for separating blood cells, and for separating cells from tissues such as bone marrow, liver, lung, spleen and solid tumors. Elutriation has many advantages: isotonic

media can be used to avoid osmotic or toxic effects of gradient materials, and there is no exposure to the effects of pelleting. Viability, purity, and recovery rates can range close to 100 percent, depending on the type of cell. And separation times are relatively brief.


The complete Elutriator System

consists of an Elutriator Rotor used in either a high speed

Beckman J-21 series centrifuge, or a low-speed J-6 series centrifuge. The new "M" versions (J2-21M and J-6M) are ideal for elutriation since they have extra capabilities built-in.

Our current applications literature describes the studies which have used the Elutriator System. Be sure to send for it if you're involved in any way with cell separations.

Write Beckman Instruments, Inc., Spinco Division, P.O. Box 10200, Palo Alto, CA 94304.

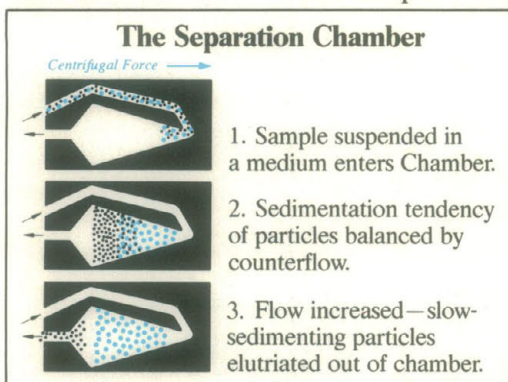


BECKMAN

Circle No. 347 on Readers' Service Card

New Translation System for Protein Processing

NEN introduces the first cell-free Translation System designed for the study of processed proteins.

Everything is in the package, including dog pancreatic microsomal membranes, dilution buffer, a positive mRNA control for monitoring system performance, and a protocol. Now you can spend less time preparing materials and free more time for evaluating results.

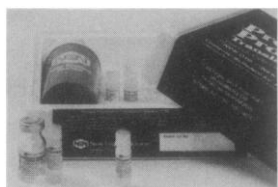

The Protein Processing Translation System, [^{35}S]- joins our family of pre-tested systems

for molecular biology in offering both convenience and reproducibility. The components, listed below, are subjected to a complete processing assay prior to shipment. The resulting Methionine, [^{35}S]- incorporation curve and autoradiogram are included in the quality control results you receive with the system.

Microsomal membranes
Microsomal membrane dilution buffer
Reticulocyte lysate
Translation cocktail
Methionine, L-[^{35}S]-
Control mRNA
Salts and distilled water

To receive the information you'll need to evaluate this exclusive new system, simply circle the inquiry number below.

Circle No. 355 on Readers' Service Card

Other Pretested NEN Systems for Biotechnology

For more information on these systems, please circle the inquiry number following the list.

Nick Translation Systems
dCTP, [^{32}P]- and dTTP, [^3H]-

3' End Labeling System
3'-dATP, [α - ^{32}P]-

5' End Labeling System
ATP, [γ - ^{32}P]-

DNA Sequencing System

Transcription Systems
UTP, [α - ^{32}P]- and GTP, [α - ^{32}P]-

Translation Systems
Methionine, L-[^{35}S]-
Leucine, L-[^3H]-
and Proline, L-[^3H]-

Circle No. 356 on Readers' Service Card

Not for use in humans or clinical diagnosis

New England Nuclear
549 Albany Street, Boston, MA 02118
Call toll free: 800-225-1572, Telex: 94-0996
Mass and Internat'l: 617-482-9595
Europe: NEN Chemicals GmbH, D-6072, W. Germany
Postfach 401240, Tel. (06103) 803-0, Telex 4-17993 NEN D
NEN Canada: 2453 46th Avenue, Lachine, Que H8T 3C9
Tel. 514-636-4971, Telex 05-821808 © 1982 NEN

NEN New England Nuclear®
a Du Pont company

LEADER IN
BIOTECHNOLOGY

LETTERS	The Missile Debate: <i>J. Bowman; D. P. Barash; M. Roeder; K. L. Lockwood;</i> Fossil Maize Pollen in Mexico: <i>P. B. Sears</i>	932
EDITORIAL	Essential Federal Information Services.....	937
ARTICLES	Solar Flares, Proton Showers, and the Space Shuttle: <i>D. M. Rust</i>	939
	Suppression of Transcription Termination by Phage Lambda: <i>D. F. Ward and M. E. Gottesman</i>	946
	Batán Grande: A Prehistoric Metallurgical Center in Peru: <i>I. Shimada, S. Epstein, A. K. Craig</i>	952
NEWS AND COMMENT	The Academic-Industrial Complex	960
	Science Education—Déjà Vu Revised	963
	Journals: Fearing the Electronic Future	964
	<i>Briefing:</i> Poland, United States Exchange Expulsions; Wyngaarden Meets the Press; More on Yellow Rain	966
	Bald Eagles on the Rise	968
RESEARCH NEWS	Is RNA Copied into DNA by Mammalian Cells?	969
	New Code Is Broken	971
	Odd Amino Acids in a Meteorite	972
	New Evidence Fuels Antarctic Ice Debate	973
ANNUAL MEETING	Call for Symposium Proposals.....	975
BOOK REVIEWS	A Treatise on the Family, <i>reviewed by W. C. Sanderson</i> ; The Biological Aspects of Rare Plant Conservation, <i>G. T. Prance</i> ; Parasitological Topics, <i>A. G. Chabaud</i> ; Ion Transport by Epithelia, <i>K. R. Spring</i> ; Theory and Experiment in Gravitational Physics, <i>D. H. Douglass</i> and <i>M. Bocko</i> ; Books Received	976
REPORTS	Hydrocarbon Resources of the Eastern Overthrust Belt: A Realistic Evaluation: <i>R. D. Hatcher, Jr.</i>	980

BOARD OF DIRECTORS

D. ALLAN BROMLEY Retiring President, Chairman	E. MARGARET BURBIDGE President	ANNA J. HARRISON President-Elect	LAWRENCE BOGORAD EDWARD E. DAVID, JR.	NANCIE L. GONZALEZ DAVID A. HAMBURG
--	-----------------------------------	-------------------------------------	--	--

CHAIRMEN AND SECRETARIES OF AAAS SECTIONS

MATHEMATICS (A) Felix E. Browder Lynn Arthur Steen	PHYSICS (B) Donald N. Langenberg Rolf M. Sinclair	CHEMISTRY (C) Charles G. Overberger William L. Jolly	ASTRONOMY (D) Irwin I. Shapiro Donat G. Wentzel
PSYCHOLOGY (J) Eleanor J. Gibson Bert F. Green	SOCIAL, ECONOMIC, AND POLITICAL SCIENCES (K) Thomas C. Schelling David L. Sills		HISTORY AND PHILOSOPHY OF SCIENCE (L) ENGINEERING (M) Erwin N. Hiebert David L. Hull Robert W. Dunlap W. Edward Lear
EDUCATION (Q) Elaine W. Ledbetter Roger G. Olstad	DENTISTRY (R) Paul Goldhaber Harold M. Fullmer	PHARMACEUTICAL SCIENCES (S) Louis A. Luzzi Robert A. Wiley	INFORMATION, COMPUTING, AND COMMUNICATION (T) Marilyn C. Bracken Madeline M. Henderson

DIVISIONS

ALASKA DIVISION		PACIFIC DIVISION		SOUTHWESTERN AND ROCKY MOUNTAIN DIVISION	
Vera Alexander President	T. Neil Davis Executive Secretary	Robert I. Bowman President	Alan E. Leviton Executive Director	Klaus D. Timmerhaus President	M. Michelle Baic Executive Officer

SCIENCE is published weekly on Friday, except the last week in December, by the American Association for the Advancement of Science, 1515 Massachusetts Avenue, NW, Washington, D.C. 20005. Second-class postage (publication No. 484460) paid at Washington, D.C., and at an additional entry. Now combined with *The Science Monthly*. Copyright © 1982 by the American Association for the Advancement of Science. Domestic individual membership and subscription (51 issues): \$48. Domestic institutional subscription (51 issues): \$85. Foreign postage extra: Canada \$24, other (surface mail) \$27, air-surface via Amsterdam \$55. First class, airmail, school-year, and student rates on request. Single copies \$2.50 (\$3 by mail); back issues \$3 (\$3.50 by mail); classroom rates on request. **Change of address:** allow 6 weeks, giving old and new addresses and seven-digit account number. **Postmaster:** Send Form 3579 to *Science*, 1515 Massachusetts Avenue, NW, Washington, D.C. 20005. *Science* is indexed in the *Reader's Guide to Periodical Literature* and in several specialized indexes.

Photoprotection by Carotenoids During Photosynthesis: Motional Dependence of Intramolecular Energy Transfer: <i>A. L. Moore et al.</i>	982
Carbyne Forms of Carbon: Do They Exist?: <i>P. P. K. Smith and P. R. Buseck</i>	984
Geochronology of Type Uquian (Late Cenozoic) Land Mammal Age, Argentina: <i>L. G. Marshall et al.</i>	986
Catastrophic Decline of a Top Carnivore in the Gulf of California Rocky Intertidal Zone: <i>M. L. Dungan, T. E. Miller, D. A. Thomson</i>	989
Bismuth Vanadate: A High-Pressure, High-Temperature Crystallographic Study of the Ferroelastic-Paraelastic Transition: <i>R. M. Hazen and J. W. E. Mariathasan</i>	991
Stem Rust of Wheat 3300 Years Old Found in Israel: <i>M. E. Kislev</i>	993
H-2 Histocompatibility Region: Influence on the Murine Glucocorticoid Receptor and Its Response: <i>C. Gupta and A. Goldman</i>	994
Human Myeloma Cell Line Carrying a Philadelphia Chromosome: <i>A. Karpas, P. Fischer, D. Swirsky</i>	997
Ribosomal Crystalline Arrays of Large Subunits from <i>Escherichia coli</i> : <i>M. W. Clark, K. Leonard, J. A. Lake</i>	999
Spectral Character of Sunlight Modulates Photosynthesis of Previtamin-D ₃ and Its Photoisomers in Human Skin: <i>J. A. MacLaughlin, R. R. Anderson, M. F. Holick</i>	1001
Decreased Nocturnal Plasma Melatonin Peak in Patients with Estrogen Receptor Positive Breast Cancer: <i>L. Tamarkin et al.</i>	1003
Insulin Stimulation of Nucleoside Triphosphatase Activity in Isolated Nuclear Envelopes: <i>F. Purrello et al.</i>	1005
Corticotropin-Releasing Factor Stimulates Accumulation of Adenosine 3',5'-Monophosphate in Rat Pituitary Corticotrophs: <i>F. Labrie et al.</i>	1007
Intracellular pH of Mitogen-Stimulated Lymphocytes: <i>D. F. Gerson, H. Kiefer, W. Eufe</i>	1009
Inability of Oxytocin to Activate Pyruvate Dehydrogenase in the Brattleboro Rat: <i>K. Hanif et al.</i>	1010
Serotonin and Octopamine in the Nematode <i>Caenorhabditis elegans</i> : <i>H. R. Horvitz et al.</i>	1012
Physiological Basis for Swimming Endurance Differences Between LDH-B Genotypes of <i>Fundulus heteroclitus</i> : <i>L. DiMichele and D. A. Powers</i>	1014
Crystallization of the Tetramer of Histones H3 and H4: <i>E. Lattman et al.</i>	1016
Epidermis of Human Skin: Pyroelectric and Piezoelectric Sensor Layer: <i>H. Athenstaedt, H. Claussen, D. Schaper</i>	1018
Interference with Dimethylhydrazine Induction of Colon Tumors in Mice by ϵ -Aminocaproic Acid: <i>J. G. Corasanti, G. H. Hobika, G. Markus</i>	1020
PRODUCTS AND MATERIALS	
Gamma Counter; Microbial Monitor; Polarizing Microscope; Nonisotopic Immunoassay; Chemistry Analyzers; Printers; Protein Processing Translation System; Literature	1024

WALTER E. MASSEY
JOHN E. SAWYER

SHEILA E. WIDNALL
HARRIET ZUCKERMAN

WILLIAM T. GOLDEN
Treasurer

WILLIAM D. CAREY
Executive Officer

GEOLOGY AND GEOGRAPHY (E)

Arthur A. Socolow
Thomas Dutro, Jr.

MEDICAL SCIENCES (N)

Steven M. Tepperman
Leah M. Lowenstein

STATISTICS (U)

Joel R. Rosenblatt
Ezra Glaser

BIOLOGICAL SCIENCES (G)

Carl Gans
Walter Chavin

AGRICULTURE (O)

Duane Acker
Coyt T. Wilson

ATMOSPHERIC AND HYDROSPHERIC

Frederic Sanders
Glenn R. Hiest

ANTHROPOLOGY (H)

John W. Bennett
Priscilla Reining

INDUSTRIAL SCIENCE (P)

Ward J. Haas
Robert L. Stern

GENERAL (X)

Daniel Alpert
S. Fred Singer

COVER

Mature bald eagle (*Haliaeetus leucocephalus*) perched on a limb at Apgar-West Glacier, Montana. These birds assemble by the hundreds each autumn to catch land-locked salmon in McDonald Creek. This year, the United States of America celebrates the 200th year of this bird as its national symbol. See page 968. [Gary Strobel, Montana State University, Bozeman, Montana 59717]

Bristol-Myers Company
is pleased to announce that
the fifth annual
Bristol-Myers Award
for Distinguished Achievement
in Cancer Research
has been awarded to
Dr. Denis Parsons Burkitt
for his clinical observation of a unique cancer
and to
Dr. Michael Anthony Epstein
for his isolation of a new virus
in Burkitt's lymphoma tissue
April 7, 1982



**New, compact Series 1000
Electronic Balances
from Sartorius...
high-priced features for
just \$495 and \$595.**

With the Series 1000, Sartorius has created two compact, portable, battery-powered electronic balances that incorporate many features found only on larger, higher-priced models.

Model 1003 (300g capacity/0.1g accuracy) is a budget conscious \$495. Model 1020 (3000g capacity/1.0g

accuracy) is just \$595. Both of these small, rugged balances feature instant electronic taring over the entire weighing range. A fast liquid crystal display reads in grams or ounces at the flick of a switch. Oversized stainless steel pans facilitate weighing large or bulky objects. And both models have built-in protection against overloading damage, dust and spills, insuring trouble-free operation.

Series 1000 Models are ideal for

use anywhere—in the lab, plant, classroom or in the field. They're extremely lightweight and can be operated on a standard nine-volt battery (or on line current with optional AC adapter.) Precision, electronic weighing on a small scale for a small price. For literature on the Series 1000, write Sartorius Balances, Brinkmann Instruments Co., Division of Sybron Corporation, Cantiague Road, Westbury, New York 11590; or call 516/334-7500.

compact
Those remarkable Δ yellow balances.


Electronic weighing on a small scale.



sartorius

For information circle reader service number 248
For demonstration circle reader service number 249

For the meeting of the world's best minds in microbiology,
you need to attend the

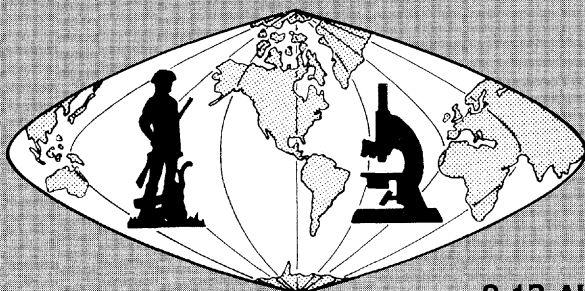

XIII INTERNATIONAL CONGRESS OF MICROBIOLOGY

International Union of Microbiological Societies

8-13 AUGUST 1982 • BOSTON, MASSACHUSETTS, USA

Sponsored by the American Society for Microbiology

The International Union of Microbiological Societies, the American Society for Microbiology, and the Organizing Committee of the Congress extend to microbiologists and other interested scientists around the world a cordial invitation to attend and participate in the upcoming XIII International Congress of Microbiology, 8-13 August, 1982. The center of Congress activities will be the John B. Hynes Veterans Auditorium in the Prudential Center, Boston. You are urged to act immediately to obtain further information. Either call the Secretary General of the Congress, at (202) 833-9680, to request additional information, or write to the address below.

Plenary Lecture:

Tumor viruses as amalgams of viral and cellular information. David Baltimore

Scientific Programs:

Over 700 contributed papers, 98 formal symposia and round tables, among four distinct programs: Interdivisional, Bacteriology, Mycology, Virology

Selected symposia include:

- Expression of plasmid genes
- Membrane transport of bacterial toxins
- Mycological research as we approach the year 2000
- Genetic control of DNA replication of phage and their hosts
- The lysogenic decision
- Beta-lactamases and the development of beta-lactam antibiotics
- Membrane energetics and transport in eucaryotic microorganisms
- Molecular basis of dimorphism
- Envelope phages as membrane models
- Current status of studies on C_1 compound metabolism
- Opportunistic infections in the compromised host
- New strategies in the control of *Neisseria* infections
- Enzyme regulation in fungi
- DNA packaging in vitro and in vivo
- Changing patterns of nosocomial infections
- Microbial fossils and traces of evolutionary changes
- Computers and microcomputers in microbiology
- Fundamental and applied aspects of ethanol production by yeasts and other microorganisms
- Cell wall macromolecules and host responses
- Structure-bound enzymes
- Fungal viruses and plasmids
- Molecular evidence of phylogenetic relationships and procaryotic evolution
- Electroimmunochemical approaches to antigenic analyses of microorganisms

Social Program:

Enjoy events such as the opening reception, symphony concert, and ballet in one of America's most historic cities; meet old friends and make new ones at this most international of meetings.

For further information, including registration package and hotel reservation forms, please write to:
Secretary General

XIII INTERNATIONAL CONGRESS OF MICROBIOLOGY

American Society for Microbiology

1913 I Street, NW • Washington, DC 20006 • U.S.A.



New from Napco®:

An autoclave engineered for the biotechnologist.

Step-saving convenience.

Completely portable unit goes anywhere. Plugs in standard 20 amp. outlet, requires no plumbing.

Timer measures only actual sterilization cycle.

New automatic precision timer starts running only after sterilization temperature is reached.

Fully automatic precision performance.

No attention is required by the operator through the complete sterilizing cycle.

Slow or fast exhaust cycle.

Switch to slow or fast exhaust cycle with the


flip of a switch.

There's much more. Send for the new Napco autoclave brochure for complete information.

napco

National Appliance Company
Subsidiary of
Heinicke Instruments Company
10855 S.W. Greenberg Road,
Portland, Oregon 97223
800-547-2555
or (503) 639-3161
Telex: 151144

Napco instruments are designed to work. But if you should ever have an equipment problem, Napco technical service assistance is available, dial toll-free 800-547-2555.

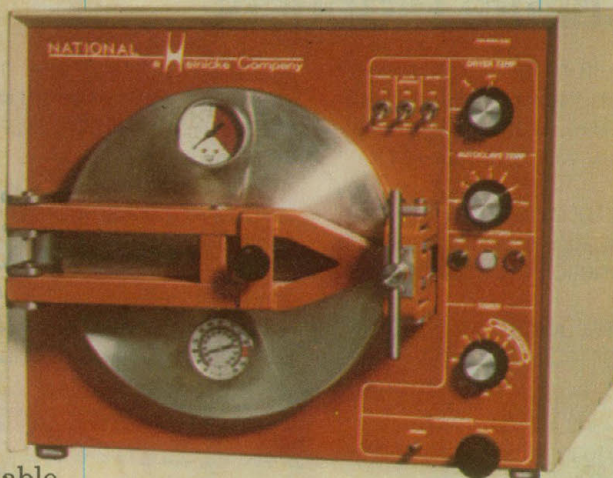

Ask your dealer's salesman today:
American Scientific Products,
Curtin-Matheson,
Fisher Scientific Co., Preiser,
Sargent-Welch, S.G.A. Scientific,
Arthur H. Thomas Co.,
In Canada:
Fisher Scientific Co., Ltd.

We put more in, you get more out.



Falcon Tissue Culture Products.

Whatever you need in tissue culture labware, Falcon probably has it. There's so much more in our disposable tissue culture products, pipets, tubes, filters and general labware because, when laboratory workers talked to us, we listened.

We applied our particular brand of innovative thinking to their problems in such specializations as hybridomas, clinical immunology/virology, genetic engineering and traditional tissue culture.

We pioneered homogeneous tissue culture treatment on plastics, and are a leader in the use of

gamma irradiation and highest quality plastics. And along the way we've grown to be the undisputed leader in an industry that accepts only the highest standards.

More new products are coming, all made with the same Falcon care, craftsmanship and rigorous quality control. So when you're looking for something in tissue culture, look at Falcon first. You may not have to look anywhere else.

Falcon Tissue Culture Products, from Becton Dickinson. Call us and let's talk.

Becton Dickinson Labware

1950 Williams Drive, Oxnard, California 93030.
Call toll-free (800) 235-5953 (outside California);
(805) 485-8711 (collect, in California).

Falcon is a registered trademark of Becton Dickinson and Company.

**BECTON
DICKINSON**

Circle No. 227 on Readers' Service Card

©1981 Becton Dickinson and Company.

New! A 4-Week Precision Delivery System.



New ALZET® 2ML4 long-duration osmotic pump makes him a super subject.

In fact, the 4-week 2ML4 will give you better results from any research animal, rat size or larger. Because our 2ML4 is the first system that delivers exactly $2.5 \mu\text{l/hr.}$, every hour, for 4 weeks.


Continuous infusion ends stress effects from repeated injections of your bioactive agent. And it's the perfect solution to variability encountered when injecting short half-life agents.

What's more, four other ALZET® osmotic pumps give you the same benefits at different delivery rates over one and two week durations.

More evidence for ALZET® effectiveness.

There are now over 400 publications in the rapidly growing ALZET® bibliography. To receive a copy, or for more information on ALZET® osmotic pumps, return the coupon below to: ALZA, Dept. Y, 950 Page Mill Road, Palo Alto, CA 94304. Or call the

technical services desk at 415-494-5323. To order, write the address above or call 800-227-9953. In California and Hawaii, call collect 415-494-5067.

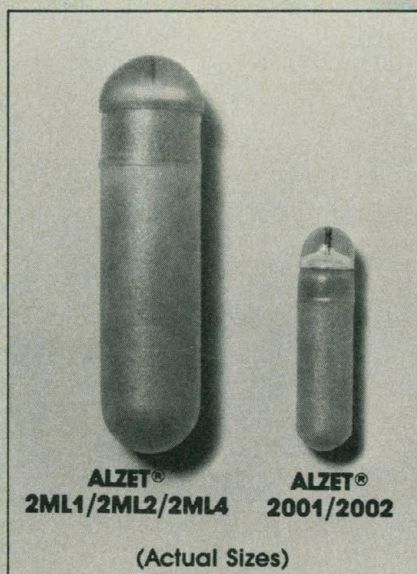

International Orders: Product inquiries, literature requests for ALZET® osmotic pumps and orders from Europe and Japan may be directed to: Scientific Marketing Assoc., 37 Mildmay Grove, London N1 4RH, England, Telephone: 01-359-5357; Tokai Irika Co., Ltd., Kurihara Building, 2-12, 3-chome, Uchikanda, Chiyodaku, Tokyo, 101, Japan, Telephone: 03-254-0052.

Circle No. 191 on Readers' Service Card

S-5

Please send me technical information on your:

☐ New ALZET® 2ML4 osmotic pump.

☐ All ALZET® osmotic pumps.

☐ Send me your new Technical Applications File and bibliography.

Name

Affiliation

Address


City

State Zip

alzet®
OSMOTIC PUMP

Computerized Cell Sorting and Analysis from Ortho

Ortho offers unexcelled systems for immunobiology researchers

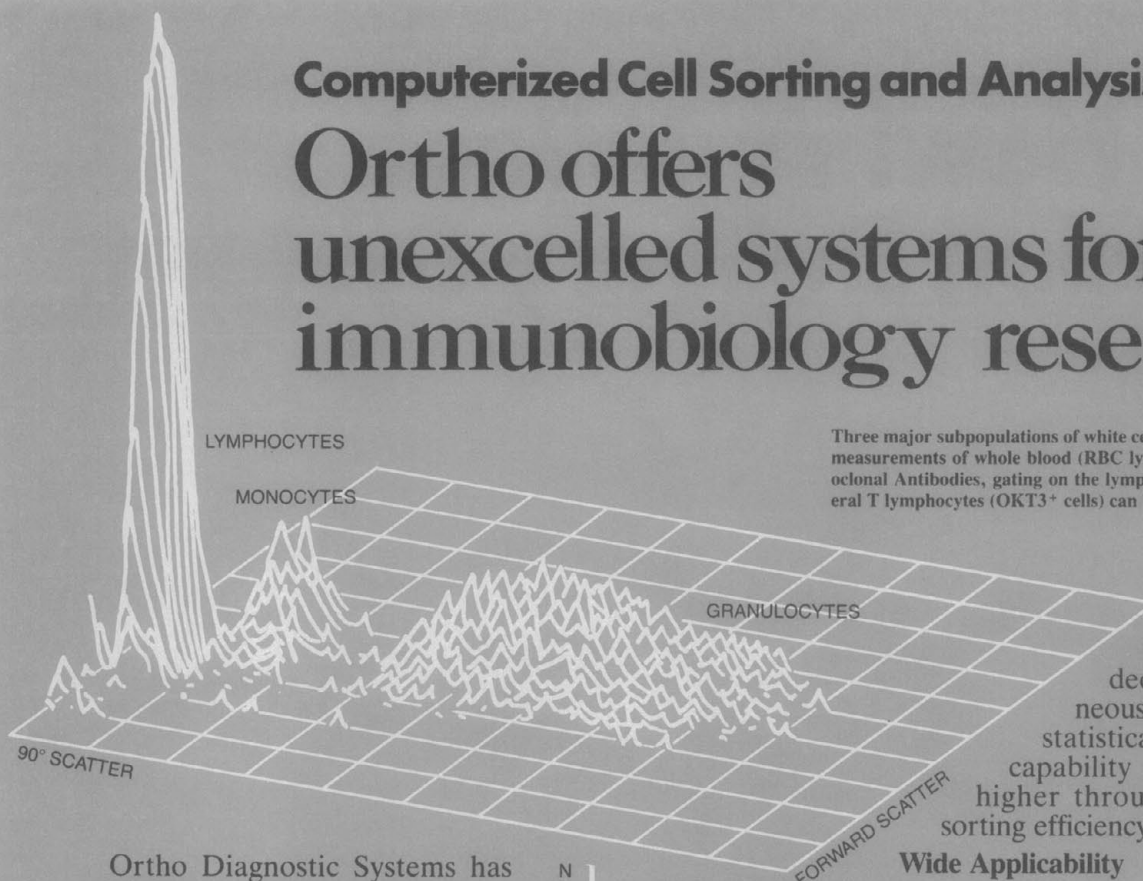

Three major subpopulations of white cells by forward and 90° light scatter measurements of whole blood (RBC lysis only). Using Ortho-mune Monoclonal Antibodies, gating on the lymphocyte population, human peripheral T lymphocytes (OKT3⁺ cells) can be quantitated and/or sorted.

it can perform rapid multiparameter sort decisions while simultaneously providing complete statistical data analysis. This capability provides significantly higher throughput and increased sorting efficiency.


Ortho Diagnostic Systems has established a firm leadership position in cell sorting and analysis instrumentation—in fact, acceptance of Ortho systems by the scientific community has been so great that these instruments now are considered the standard of performance in the field.

Versatile Instrumentation

The Ortho Cytofluorograf* Flow Cytometer is an example of that leadership and commitment to research. As the only flow cytometer to offer a choice of five lasers, (including two installed as standard), the Cytofluorograf provides exceptional applications flexibility.

Flexibility is enhanced further by a modular design which allows the system to be upgraded as research needs change. Performance too is outstanding with analysis rates of up to 10,000 cells per second. And only Cytofluorograf offers direct measurement of seven optical parameters including axial light loss, a valuable technique that allows the use of non-fluorescent stains.

The new 2150 Computer and software can acquire up to 16 histograms and simultaneously display any four of these, all in real-time, allowing the operator complete control over all aspects of the experiment. Because 2150 has two independent microprocessors,

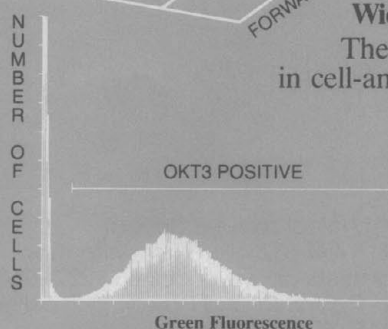

Wide Applicability

The system has proved extremely valuable in cell-antigen and tumor studies. It can also be employed for research into cell sizing, DNA and chromosome analysis and many aspects of cell kinetics and microbiology. The exceptional versatility of the system allows it to keep pace with the researcher's interests, never to limit them.

Total Professional Support

Ortho's interest in cellular biology extends beyond producing instrumentation. For example, Ortho has developed Ortho-mune* monoclonal antibodies, the first panel (for research use only) of murine monoclonal antibodies reactive with human T lymphocytes. The company maintains well-equipped applications and training laboratories to develop new techniques and to provide consultative services for customers.

The results of these investigations and others are published in Protocols so that other researchers may share the results. Free copies are available by contacting Ortho.

Ask about the Cytofluorograf Flow Cytometer, 2150 and the complete range of Ortho services. Contact Ortho Diagnostic Systems Inc., Research Instruments Division, 410 University Avenue, Westwood, MA 02090 (617) 329-6100.




Ortho Diagnostic Systems Inc.

* T.M. of Ortho Diagnostic Systems, Inc.

Copyright 1981 Ortho Diagnostic Systems, Inc.

Circle No. 214 on Readers' Service Card

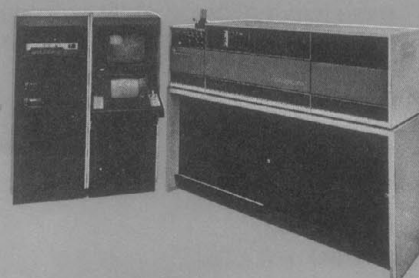

ANNOUNCING: **a series of COULTER CLONE™ Antibodies**

The most comprehensive panel of monoclonals for investigating such conditions as:

- Immune deficiencies
- Autoimmune disorders
- Non-Hodgkins leukemia and lymphomas
- Monocytic and granulocytic leukemia
- Lymphocytic and lymphoblastic leukemia
- T & B lymphocyte subclassifications
- Suppressor/helper ratios

Unprecedented specificity for defining certain disease states.

Highest labeling density available.

Complete functional measurements.

Unique CALLA, I2 (HLA - D/DR) and B1 cell markers.

SPECIAL INTRODUCTORY OFFER

We are making seven of our monoclonal antibodies available now: T11, T4, T8, B1, I2, J5, and Mo2. For this offer, they are being specially packaged in 25-test vials at a price of only \$37.50 per vial. The offer is limited to one vial per customer of any or each type.

Call us toll free at 800-327-6531, extension 501 for more information, extension 359 to order. Or you can complete and mail the coupon to:

S

name and title

company/lab


street address

city

state, zip


phone

Circle the antibodies of interest to you.
T11 T4 T8 B1 I2 J5 (CALLA) Mo2

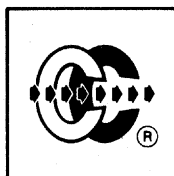

Marketing Department
COULTER ELECTRONICS, INC. 590 West 20th Street Hialeah, FL 33010
Circle No. 226 on Readers' Service Card

For research use only, not for diagnosis

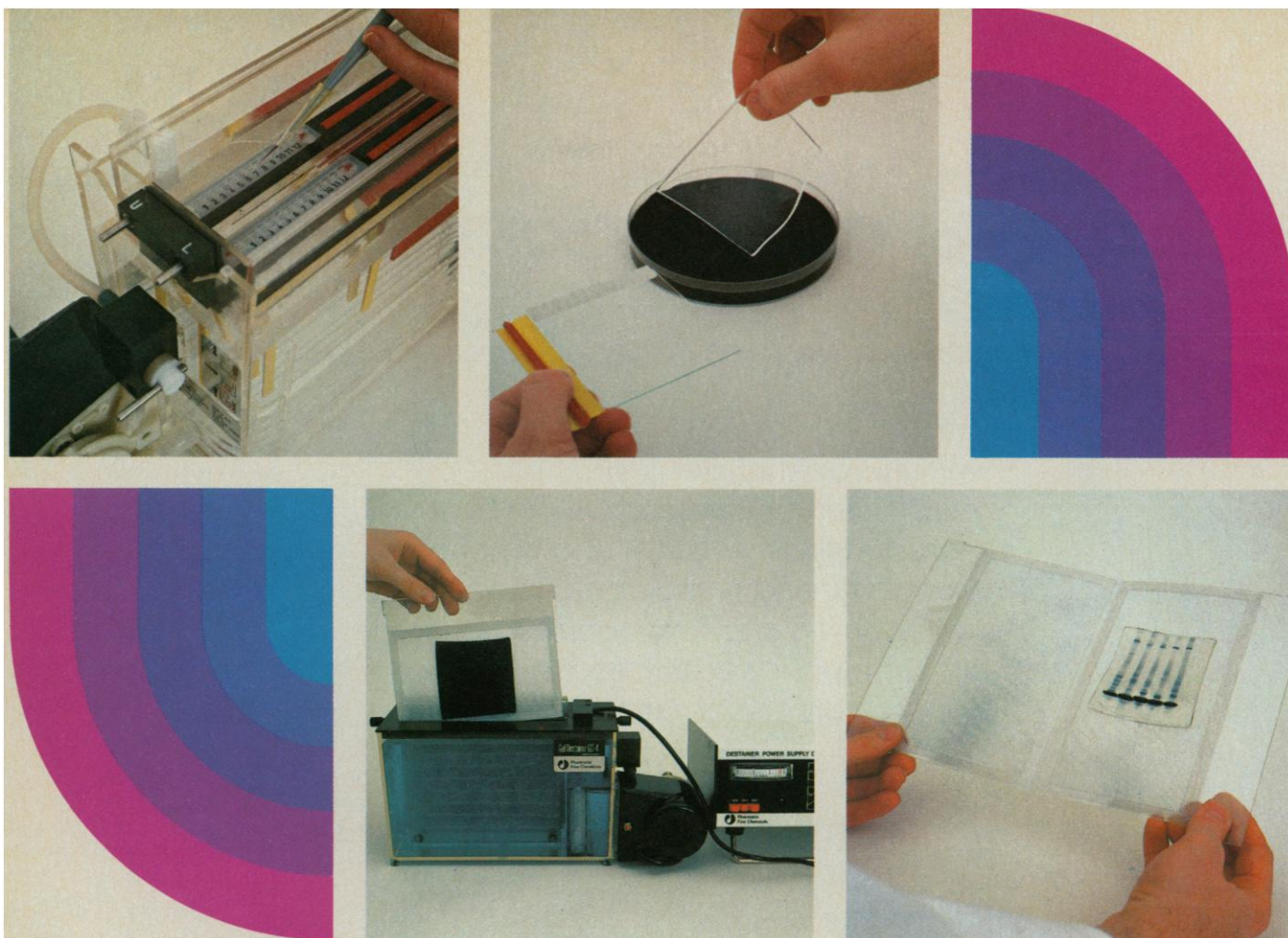


2D or Not 2D?


The choice is yours... using only one electrophoretic unit.

Important gel electrophoresis methods, such as those established by O'Farrell, Weber and Osborn, Anderson and Anderson, and Laemmli, require many gel configurations and conditions. Gel rods and slabs of varying thickness, width, and length are essential to perform all these techniques.

The GE-2/4 and GE-2/4 LS are electrophoretic units designed for a wide variety of techniques including two-dimensional work, gradient gel electrophoresis, SDS-PAGE, and isoelectric focusing in gel rods. Each unit is designed to accept a wide variety of sizes of either gel rods or gel slabs.

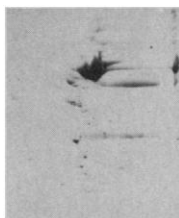

A 2-D gel (Anderson and Anderson) of human serum run in a GE-2/4 LS using Pharmalyte® pH 3-10 in the first dimension and a 10-20% gradient gel in the second dimension.

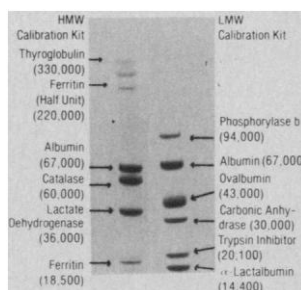

Proteins run on a Pharmacia PAA 4/30 gradient gel in 0.2% SDS.


Send for more information on the system that gives you YOUR CHOICE.

Pharmacia Fine Chemicals
Division of Pharmacia, Inc.
Piscataway, New Jersey 08854
Phone (201) 469-1222

See us at the TCA Meeting in San Diego—Booth #6

Pharmacia
Fine Chemicals

Drummond Liquid Handling Systems ...

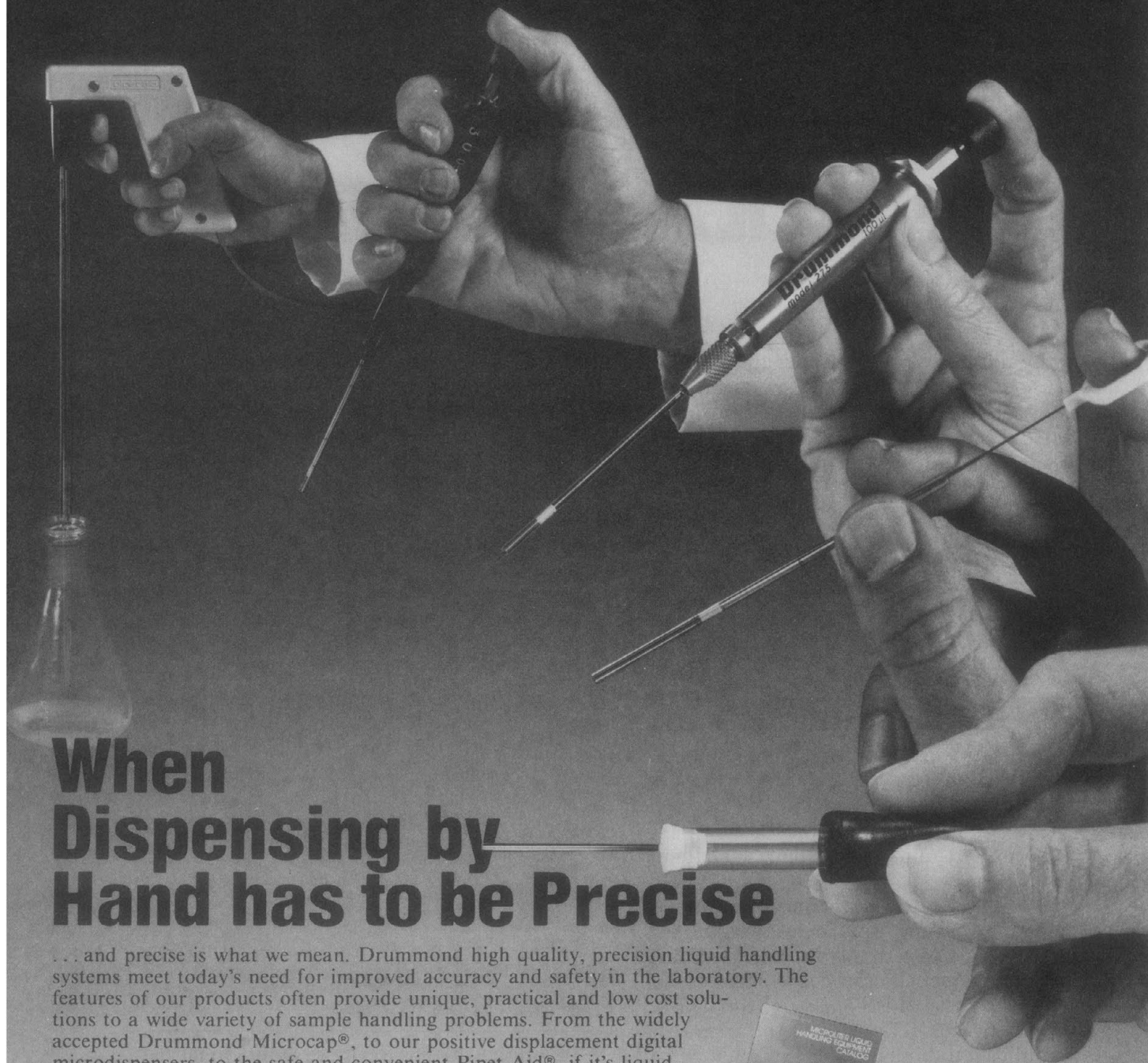

When Dispensing by Hand has to be Precise


... and precise is what we mean. Drummond high quality, precision liquid handling systems meet today's need for improved accuracy and safety in the laboratory. The features of our products often provide unique, practical and low cost solutions to a wide variety of sample handling problems. From the widely accepted Drummond Microcap®, to our positive displacement digital microdispensers, to the safe and convenient Pipet-Aid®, if it's liquid, there's a Drummond product that can handle it ... better.



DRUMMOND
SCIENTIFIC COMPANY

500 PARKWAY, BOX D □ BROOMALL, PA 19008
Call toll-free: 800/523-7480 or in Pennsylvania call collect 215/353-0200

Circle No. 234 on Readers' Service Card

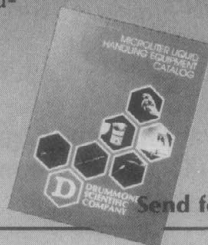

Send for Catalog

Ortho-mune^{*} monoclonal antibodies in only two years generated more than 200 published papers

In a way, this is a tribute to Köhler and Milstein. Because they're the ones who gave us the idea. But the scientists and engineers at Ortho deserve credit, too, for making it happen.

It started, of course, with monoclonal antibodies. Six years ago, when George Köhler and Cesar Milstein showed the world how hybridomas could be made to produce antibodies against specific immunogens, there was no way to predict exactly which questions the new technology would answer.

OKT^{*}: WELL-CHARACTERIZED MONOCLONAL ANTIBODIES[†]

The questions we had at Ortho were immunological ones. Our researchers were particularly interested in human lymphocytes. We developed a series of monoclonal antibodies directed against purified human T cells and called them the OKT panel.

Then the real work began. Because clear-cut results can come only from well-characterized reagents, scientists at Ortho put a lot of emphasis on defining the specificity of each clone. Input also came from scientists outside Ortho who worked with the OKT reagents. And as more has been published about the reagents, they've become more and more useful.

New potential applications came to mind. It began to look like the OKT reagents could

standardize T lymphocyte subset analysis. Some immune deficiency disease states turned out to be associated with abnormal ratios of T cell subsets. Research oncologists have now begun using the OKT reagents to characterize proliferative lymphocytic disease states, and they have been used to study immune responses to renal transplantation.




MONOCYTES AND Ia-POSITIVE CELLS

We continued to develop more monoclonal antibodies. The new ones could differentiate monocytes and Ia-positive cell subsets. More applications were developed—sorting human peripheral blood cells, elimination of cell populations by complement lysis, studies on cytotoxicity and histochemical staining of tissue sections.

^{*}Trademark

[†]For research use only.
Not to be used for diagnostic procedures.

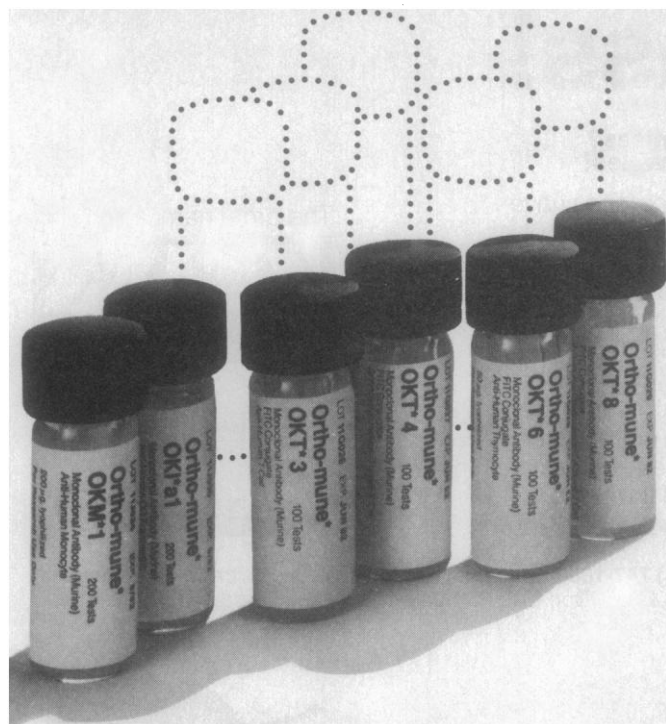

FITC-CONJUGATED ANTI-HUMAN T CELL REAGENTS

We know how busy you are. Our new fluorescein-conjugated OKT monoclonal antibodies should help. They reduce staining procedures to one time-saving step. And more new monoclonal antibodies and fluorescent reagents are being developed and characterized at Ortho.

THE LEADING EDGE

Applications for the Ortho-mune products are following quickly on the heels of research developments. Research oncologists, rheumatologists and immunologists in and out of Ortho are using the System to take a closer look at leukemia, rheumatoid arthritis, multiple sclerosis, myasthenia gravis and response to infections. We've become more than just producers of monoclonal antibodies. We have a commitment to development of integrated systems for analysis of the human immune system. Ortho Immunology Systems is a name that seems to cover it.

Ortho brings you innovative technology—reagents and instruments developed by research immunologists for immunology research.

We'll be pleased to talk with you about setting up a system for lymphocyte subset analysis in your own laboratory. Call us at (800) 526-3875. In New Jersey, call (201) 524-1365.

CELL SORTING AND ANALYSIS: CYTOFLUOROGRAF®


We created the Cytofluorograf to meet the growing demands for laser based cell analysis and sorting. Then we combined it with the Model 2150 computer system to provide the flexibility needed for ever-changing research demands in immunobiology, chromosome analysis and molecular genetics. That's the kind of innovation that makes the Cytofluorograf a leader in flow cytometry.

ORTHO SPECTRUM III*

If you are involved in large volume clinical research projects, you may need high speed analysis of whole blood, without the inconvenience of time-consuming sample preparation. That's what the Ortho Spectrum III*. Integrating Ortho-mune® monoclonal antibody technology and laser flow cytometry, this flexible microcomputer-based system can be programmed for accurate cell differentiation on the basis of size and immunofluorescence.



Ortho Diagnostic Systems Inc.

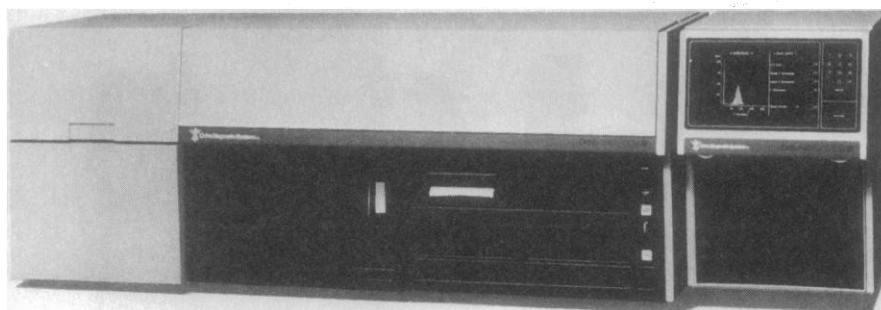

*Trademark

Formerly Orthoclone

Head of the class!

The Zeiss Universal does more

When it comes to "head-on" competition, the others can't touch Zeiss Universal Microscopes. Built into the tube head of the Universal and UEM are the optional "Optovar" magnification changer; a completely flexible beamsplitter that will send all the light to the observer, or all to a camera, or divide it; and slots for analyzers, filters, wave plates and compensators.

Promises...promises? Deliveries...deliveries!


With Universals you won't have to wait for the development of accessories. Everything needed for every technique in transmitted and reflected light already exists and is proven in use. And, of course, the quality of Zeiss optics has been the goal toward which others have been striving for generations. Today, Zeiss optics are still at the head of the class.

Quality service—Expert dealers.

The great name in optics



Carl Zeiss, Inc., One Thornwood Drive, Thornwood, N.Y. 10594 (914) 747-1800. Branches: Atlanta, Boston, Chicago, Houston, Los Angeles, San Francisco, Washington, D.C. In Canada: 45 Valleybrook Drive, Don Mills, Ontario, M3B, 2S6. Or call (416) 449-4660.

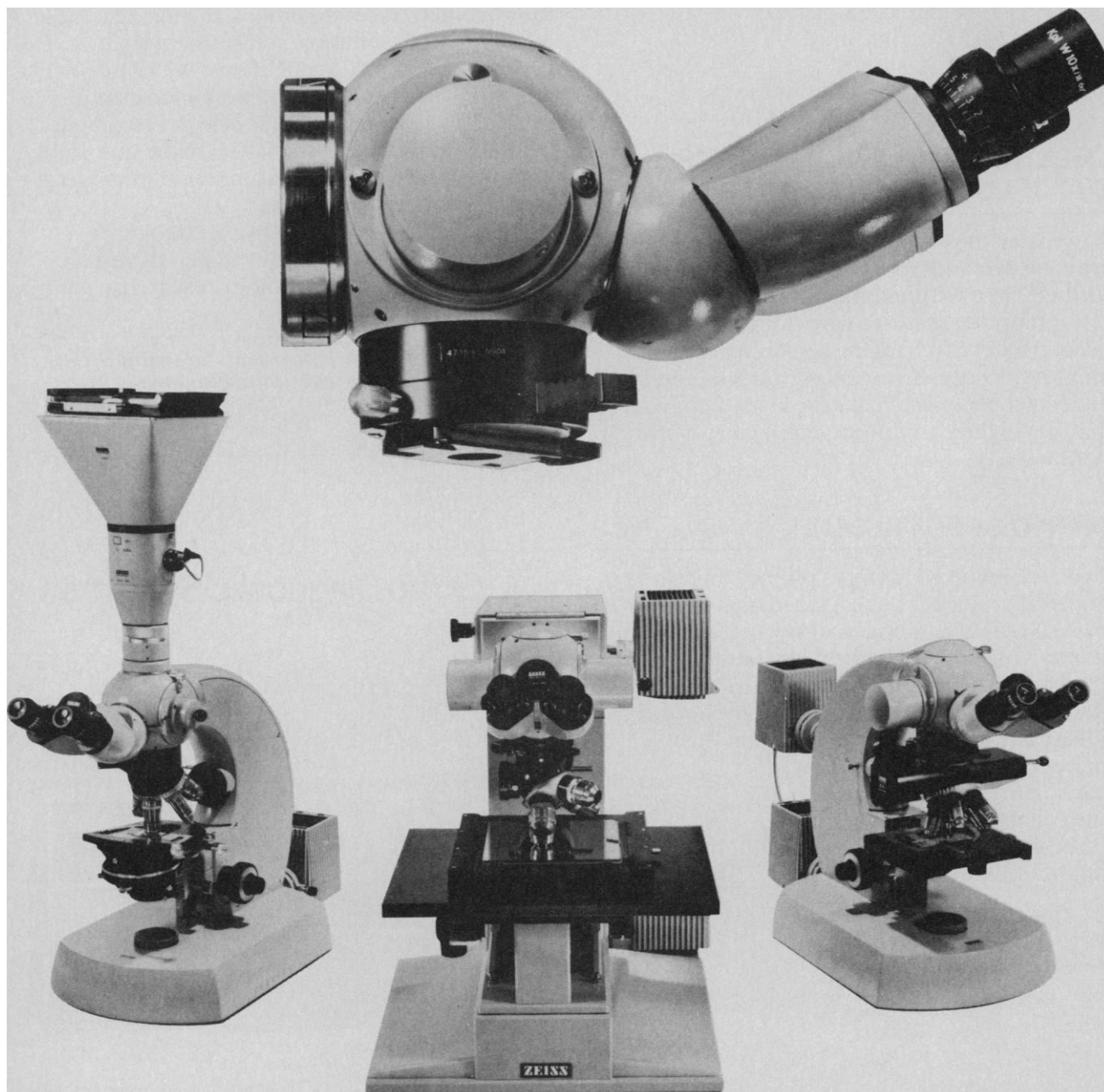

FOR TODAY'S FACULTY AND COLLEGE STAFF MEMBERS* FROM 18 TO 80.

Whether you're thinking retirement or not, review the plan that provides for cash withdrawal and/or a lifetime income.

TIAA-CREF Supplemental Retirement Annuities (SRAs) offer you substantial flexibility including cash withdrawal and/or lifetime retirement income. You can even reduce your income taxes now!

You can begin contributions to an SRA at any age and begin benefits at any age up to age 71 unless you are still employed (then you can delay beginning benefits until age 80). For example, you could start contributions at age 25, and choose to begin benefits or withdraw cash at age 34, 40 or 50, regardless of your employment status.

Get your money at any time.

You can receive benefits as a lifetime income or over a fixed period of from 2 to 10 years. What's more, if you need it (even while employed by your current employer), you can withdraw all the money you have accumulated by surrendering your contracts. Or, you can withdraw \$1,000 or more every six months. There is never a cash surrender charge.

Contributions are tax-deferred, so you pay less income taxes now.

The federal income tax on your contributions is deferred until they are paid to you as benefits. So, you pay less tax now.

Changing employers? Take SRAs with you.


Since you own your Supplemental Retirement Annuities, you take them with you if you leave your current employer. You can make contributions through any institution that makes Supplemental Retirement Annuities available to staff members. Contributions can be as little as \$25 a month.

Full information.

Complete and mail the coupon for an SRA Information Kit today. You'll get full details about all the advantages SRAs have to offer, why this plan suits so many financial situations and age groups and how much you may contribute to the plan.

*TIAA-CREF provides annuities and other services for employees of colleges, universities, private schools and certain other nonprofit tax-exempt educational and research institutions.

HELP YOURSELF TO A BRIGHTER FINANCIAL FUTURE ... SEND FOR A FREE INFORMATION KIT.

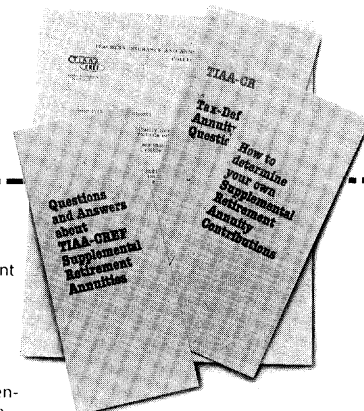

Teachers Insurance and
Annuity Association of
America—College Retirement
Equities Fund
730 Third Avenue
New York, N.Y. 10017

Please send me full details
about TIAA-CREF Supplemental
Retirement Annuities, the
flexible tax-deferred annuity plan that offers the opportunity to accumu-
late funds for additional retirement income and the option for cash
withdrawal.

Name _____

Address _____

City _____


State _____ Zip _____ Date of Birth _____

Name of Institution _____

Please let us know if you are participating in a TIAA-CREF retirement plan
at your institution. ☐ YES ☐ NO

S52882

How to get a headstart on new projects.

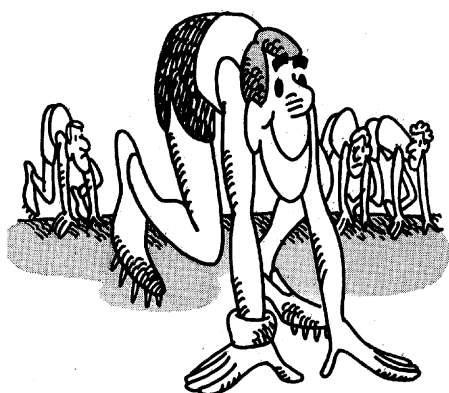

There's a lot of information scattered around the world that can help you, if only you can find it. That's the job for Lockheed's Dialog, the world's leading online information retrieval system.

This computerized sleuth puts you in instant touch with more than 130 databases and provides abstracts of journal articles, conference papers, research reports, patents and news stories—and all of this information goes back 15 years.

Dialog covers virtually every field, from acoustics to zero gravity.

All you need to use Dialog is a standard computer terminal and a telephone. Your library probably subscribes to Dialog.

Involved in some phase of a new project? Better get involved first with Dialog. It can save you time, money and effort—and lots of each.

For more information, contact Dialog Information Services, Dept. 83, 3460 Hillview Ave., Palo Alto, CA 94304. Call toll-free (800) 227-1927. In California, call (800) 982-5838.



Circle No. 358 on Readers' Service Card

LETTERS

The Missile Debate

R. Jeffrey Smith's articles on U.S. strategic thought and the MX missile (News and Comment, 2 Apr., p. 30; 9 Apr., p. 150; 16 Apr., p. 270; 23 Apr., p. 388; 30 Apr., p. 492; 7 May, p. 596; 21 May, p. 828) make a commendable contribution to the current debate over the nuclear arms freeze. These articles expose the faulty thinking behind the Reagan Administration's claims that the Soviet Union has achieved superiority over the United States in nuclear weapons.

However, Smith dismisses too lightly the question of why the Administration is trying so hard to sell the United States on a nuclear arms race on the basis of such a deception. . . . In the book *Protest and Survive* (Monthly Review Press, New York, 1981), Daniel Ellsberg argues that the Reagan and Carter administrations have been deliberately seeking a first-strike capability in order to threaten the Soviet Union with annihilation in case of a confrontation over a Third World country like Afghanistan, Iran, or El Salvador. . . . Striving for a first-strike capability today is not just a ticket for an incredibly expensive arms race in the face of the law of diminishing returns. If the Administration ever believes the first-strike capability has been attained, then the slightest miscalculation during a Third World crisis could turn into a nuclear war. That is one reason why millions of people here and in Europe are taking actions to stop the nuclear arms race.

JOSEPH BOWMAN

5809 Kinoll Avenue,
Cincinnati, Ohio 45213

Science and R. Jeffrey Smith are to be congratulated for the superb News and Comment series on basing modes for the MX missile. However, beyond the exorbitant expense, the environmental danger and the Strangelovian nature of the whole affair, there is another, more basic problem that is not mentioned: whereas the basing mode has received the bulk of political, military, and technical scrutiny, the MX missile itself has remained relatively unscrutinized. Ironically, the arguments against the missile—whatever its basing mode—are even more cogent than those already raised against the various plans for its deployment.

Like the planned Trident II, the MX possesses such hard-target lethality as to make it a credible first-strike weapon

and a threat to Soviet land-based intercontinental ballistic missiles (ICBM's). . . . An important commandment for survival in the nuclear age should read "Thou shalt not threaten thy neighbor's deterrent." Whether the MX is based in a vulnerable or invulnerable manner, it violates that commandment and must be opposed.

DAVID P. BARASH

Departments of Psychology and
Zoology, University of Washington,
Seattle 98095

The recent series of articles by R. Jeffrey Smith are among the clearest and most detailed analyses of the topic it has been my pleasure to read.

While reasonable persons may differ in their view of the necessity or lack thereof of an arms buildup, no thinking person should arrive at a conclusion without at least reading these articles. . . .

MARTIN ROEDER

Department of Biological Science,
Florida State University,
Tallahassee 32306

In response to Smith's article "Carter's plan for MX lives on" (News and Comment, 30 Apr., p. 492), could it be that "BMD" stands for "Bury Me Deep," where the "Me" could mean missile, me, the citizenry, or even civilization?

Or, maybe it means simply "Big Money Defense."

KARL L. LOCKWOOD

Department of Mechanical
Engineering, Ohio State University,
Columbus 43210

Fossil Maize Pollen in Mexico

The report of maize pollen 70 meters under Mexico City in a layer of Pleistocene age has been publicized (1) and discussed (2-3) enough to warrant reexamination of the evidence. At the request of its finder the identity of the pollen has been checked and confirmed (2). To suggestions that the supposed Pleistocene pollen is a contamination brought down from younger levels in coring, it can be pointed out that the pollen was found in samples taken from the intact interior of precision cores.

Overlooked meanwhile have been (i) the basis for assuming that the 70-meter level is of Pleistocene age; (ii) the effect of the well-known subsidence of the area cored; and (iii) ample evidence elsewhere in the Basin of Mexico of two

COMPUTER AIDED CHEMISTRY

The concept that is revolutionizing laboratory productivity.

Perkin-Elmer is leading a revolution in productivity, in laboratory management, in manpower use and economics...a revolution called Computer Aided Chemistry.



or part time as a computer and part time as an instrument controller. Result: A very flexible and easy-to-use small laboratory data management system that can be configured according to

LEVEL 1—INSTRUMENTS THAT PUT PEOPLE BACK TO WORK ON CHEMISTRY.

At the heart of Computer Aided Chemistry is the same tiny microprocessor that is transforming commercial and industrial electronics. Built into analytical instruments, it makes possible an awesome jump in their efficiency. Setup and calibration, methods, autosampling, sample identification and data reduction are all controlled electronically.

The broadest line. Perkin-Elmer offers the world's broadest line of microprocessor controlled analytical instruments. *Almost every analytical instrument Perkin-Elmer offers is microprocessor controlled.* If your laboratory uses one or more of them, you are already working with Computer Aided Chemistry.

LEVEL 2—DATA STATIONS MAKE INSTRUMENTS AND ANALYSTS MORE PRODUCTIVE.

The microprocessor also opens ways for instruments to interact with advanced control devices, such as Perkin-Elmer Data Stations.

A Data Station is actually a microcomputer system through which operator and instrument communicate. Programmable in BASIC, it can order and compare analyses, record and manipulate data. Results can be displayed on a CRT or as hard copy on a printer.

Result: More savings in time, another increase in productivity.

Consider the Perkin-Elmer DMS/200 Data Management System, for instance. Basic hardware is a Model 3600 Data Station. But Perkin-Elmer software enables the Data Station to function either full time as a stand-alone computer,

your needs, generate reports, search the data base, and draw off data for use in programs.

LEVEL 3—LIMS/2000 UNLOCKS THE FULL POTENTIAL OF COMPUTERIZATION.

LIMS/2000—Laboratory Information Management System—consists of advanced software and a broad array of instruments and Data Stations interacting with a powerful 32-bit Perkin-Elmer minicomputer.

LIMS/2000 not only fosters a quantum leap in productivity, it gives the manager of even the largest laboratory total control over the flow of work.

A single source. Perkin-Elmer Computer Aided Chemistry integrates instruments, computers, and software through a single, reliable source. This means you can begin a system with instruments, add Data Stations as workload grows, and a mini-computer still later *without* obsoleting any hardware or software purchased in the process.

Let Perkin-Elmer show you.

Perkin-Elmer will make a free survey of your laboratory and show you how Computer Aided Chemistry can improve productivity and efficiency. We also have a brochure on Computer Aided Chemistry.

Ask for a copy. Just write or call:



Software—The Perkin-Elmer advantage.

The catalyst for communication at every level of Computer Aided Chemistry is software. And Perkin-Elmer offers the broadest repertoire of proven ready-to-use software in the analytical field. More than 50 software programs for Computer Aided Chemistry are available now, and cover the full range of analytical disciplines.

Perkin-Elmer Corp., Analytical Instruments, Main Ave. (M.S.12), Norwalk, CT 06856 U.S.A. Tel: (203) 762-1000.

Bodenseewerk Perkin-Elmer & Co., GmbH, Postfach 1120, 7770 Ueberlingen, Federal Republic of Germany. Tel: (07551) 811.

Perkin-Elmer Ltd., Post Office Lane, Beaconsfield, Bucks HP9 1QA, England. Tel: Beaconsfield (049 46) 6161.

PERKIN-ELMER
Responsive Technology

R&D

in the

FY 1983

Budget: Impact and Challenge

7th Annual

**AAAS Colloquium on
R&D Policy**

23–24 June 1982

**The Shoreham Hotel
Washington, D.C.**

Issues include:

R&D policy and the FY 1983 budget □ *the climate for industrial R&D* □ *R&D programs of key federal agencies* □ *impacts on the scientific and engineering communities*

Participants include leaders from:

OMB, OSTP, and federal agencies □ *Congress* □ *industry* □ *the scientific and engineering communities*

In addition:

Research and Development, AAAS Report VII, by Willis H. Shapley, Albert H. Teich, and Jill P. Weinberg, will be provided in advance to Colloquium registrants. The Report covers R&D in the federal budget for FY 1983 and other topics on R&D and public policy. Registrants will also receive the published Proceedings of the conference.

For further details, write:

R&D Colloquium

**AAAS Office of Public Sector
Programs**

**1776 Massachusetts Ave., NW
Washington, D.C. 20036**



**American Association for the
Advancement of Science**

cultural horizons, Archaic 2000± to 0± B.C.) and Nahua (A.D. 800± to 1519) separated by lava and volcanic ash (4).

In 1948 and 1949, two cores, Madero and Bellas Artes, from the center of the colonial city of Mexico were made available for study at Oberlin College in the hope that they might furnish information as to conditions outside the limits of continental glaciation. The cores, 3 to 5 inches in diameter and 70± meters in depth, were generously furnished by the distinguished engineer Leonardo Zeevaert, whose studies preceded the successful construction of the 43-story Latino Americana at a site notable for subsidence (5). Kathryn Clisby analyzed the pollen (6), and at her wise insistence Fred Foreman examined the sediments (7).

So far as I can determine, belief in the Pleistocene record of the cores has been shared by all concerned and rests upon a suggestion by Kirk Bryan (8), whose caution has not been sufficiently appreciated. Since upland soils whose fossil content proves their Pleistocene origin are well known in the Basin of Mexico, it seemed reasonable to suggest that cognate horizons had accumulated as sediment in a central deep. Thus for the equivalent of a fossil soil of known Tazewell-Cary age, Bryan wrote: "*Se supone (sic) que existe en el centro de la Cuenca,*" but further, "*No se pretenden resultados absolutamente conclusivos,*" followed by a clear statement of the need for complete study.

Actually Foreman (7) found no evidence that there had ever been a deep lake; instead the presence of oölites, sand, root-holes and roots indicated a succession of swamps and shallow lakes. Puzzled, yet apparently believing in the antiquity of the cores, Foreman suggested that excess water must have drained away into permeable lake margins as a sequence of surfaces piled up. I am unable to find mention of the likelihood that the sinking of successive surfaces rather than a building up from the bottom of an original basin was involved.

Yet Zeevaert's profile (5) of seven cores shows layers sloping downward some 20 meters in 4 kilometers to their position in the two deep cores studied. Projected inward from the numerous sites (4) where the Archaic, ash or lava, and Nahua horizons occur within a few meters of the surface and given a mean annual rate of less than 4 centimeters of subsidence, the Archaic horizon could easily have sunk to 70 meters in the 2000± years since it was covered by the eruption of Xitle. The absence of arti-

facts along with the maize pollen is in no way remarkable; exploration of the 70-meter level has been limited to two slender cores; nor were the one-time swamp and shallow lake surfaces favorable to heavy accumulation of artifacts.

The pollen profiles (7) noted by Foreman as an essential complement to the sediment record show a threefold pattern. The top 6 to 8 meters record a complex of composites, amaranth (doubtless the crop *quelite*), and grasses, including maize. Along with an abundance of artifacts in the surface *relleno* or fill, this complex is a valid index of Nahua lake-margin or *chinampa* horticulture. Below it are some 50 or more meters of the old surfaces of shallow lakes and swamps mixed with weathered and unweathered volcanic debris until the 55.6-meter level is reached in Madero and the 71.2-meter level in Bellas Artes. Here the complex reappears, including maize pollen in Bellas Artes, thus reconciling the core record with that elsewhere in the basin. Nor is the recent evidence of identity of maize and tetraploid teosinte pollen (3) necessary to explain Pleistocene occurrence of apparent maize pollen because there is good evidence of the cultivation of maize before 2000 B.C. (9).

Subject always to further research it is my judgment that the pollen at 70 meters is an index of Archaic horticulture and not of wild Pleistocene maize. I deeply regret any embarrassment this may cause others but am less concerned that it makes nonsense of my efforts (10, p. 526) to relate the two Mexican profiles to Pleistocene events than satisfied that it fits the Archaic-Teotihuacan-Nahua sequence to which my esteemed friend and *tocayo* Pablo Martinez del Rio invited my attention in 1948.

PAUL B. SEARS

Las Milpas,

Taos, New Mexico 87571

References

1. M. Elting and M. Folsom, *The Mysterious Grain* (Evans-Lippincott, New York, 1967), p. 78; D. Wilson, Ed., *The New Archeology* (Meridian Reprint, New York, 1976), p. 175.
2. P. C. Mangelsdorf, E. Barghoorn, U. Banerjee, *Bot. Mus. Leaflet. Harv. Univ.* 26, 237 (1978).
3. G. W. Beadle, *Science* 213, 890 (1981); J. C. Cámara-Hernández and P. C. Mangelsdorf, *Bussey Inst. Publ.* 10 (1981), p. 3; P. C. Mangelsdorf, L. M. Roberts, J. S. Rogers, *ibid.*, p. 39.
4. P. B. Sears, *Bull. Geol. Soc. Am.* 63, 241 (1952).
5. L. Zeevaert, *Mem. Congr. Cient. Mex.* 5, 58 (1953).
6. K. Clisby and P. B. Sears, *Bull. Geol. Soc. Am.* 66, 511 (1955).
7. F. Foreman, *ibid.*, p. 475.
8. K. Bryan, *Bol. Soc. Geol. Mex.* 13, 1 (1948).
9. P. C. Mangelsdorf, R. S. MacNeish, G. R. Willey, in *Handbook of Middle American Indians*, R. Wauchope, Ed. (Univ. of Texas Press, Austin, 1965), vol. 1, pp. 427–445.
10. P. B. Sears and K. Clisby, *Bull. Geol. Soc. Am.* 66, 521 (1955).

Napco® promises you:

No other water jacketed CO₂ incubator can match the new 6000 series.

Over 15 advanced features are making incubator history. Here are a few of them.

More accurate control of chamber temperature.

By positioning sensors not only in the water jacket, *but in each separate chamber*, Napco permits more precise chamber monitoring than other water jacketed incubators.

The only water jacketed CO₂ incubator you can get really clean.



Shelf standards and plenums can be easily removed for cleaning in five minutes. This leaves an empty, perfectly smooth chamber that can be swabbed for easy management of spores and contaminating

organisms.


Unique wet-jet air flow system.

The exclusive wet-jet air flow system provides maximum chamber uniformity and assures the highest relative humidity available in any incubator.

And the list goes on. Feature after feature are designed for your convenience and precision performance. Send for the complete brochure.

All Napco water-jacketed

incubators are available in three space-saving models.

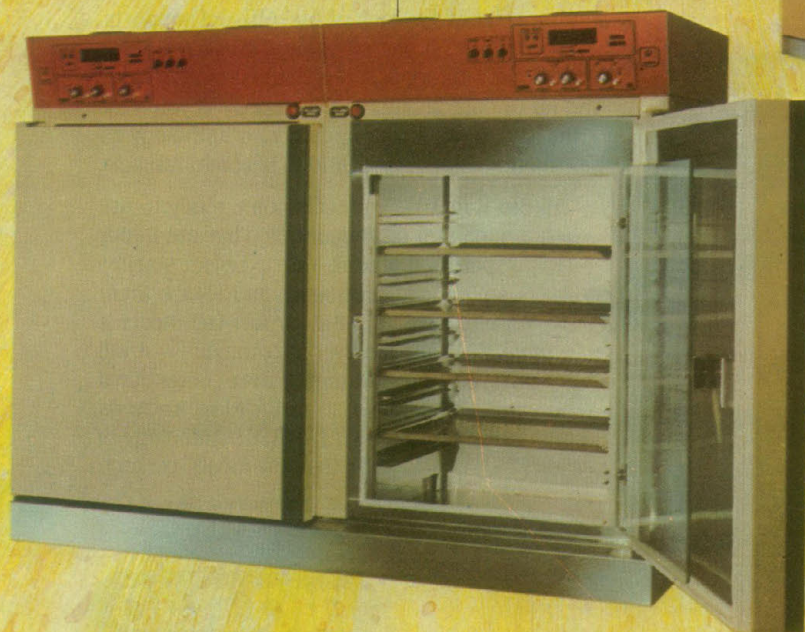


napco

National Appliance Company
Subsidiary of
Heinicke Instruments Company
10855 S.W. Greenberg Road,
Portland, Oregon 97223
800-547-2555 or (503) 639-3161
Telex: 151144

Ask your dealer's salesman today:
American Scientific Products,
Curtin-Matheson, Fisher Scientific Co.,
Preiser, Sargent-Welch, S.G.A. Scientific,
Arthur H. Thomas Co.,
In Canada: Fisher Scientific Co., Ltd.

Napco instruments are designed to work. But if you should ever have an equipment problem, Napco technical service assistance is available, dial toll-free 800-547-2555.

Circle No. 32 on Readers' Service Card

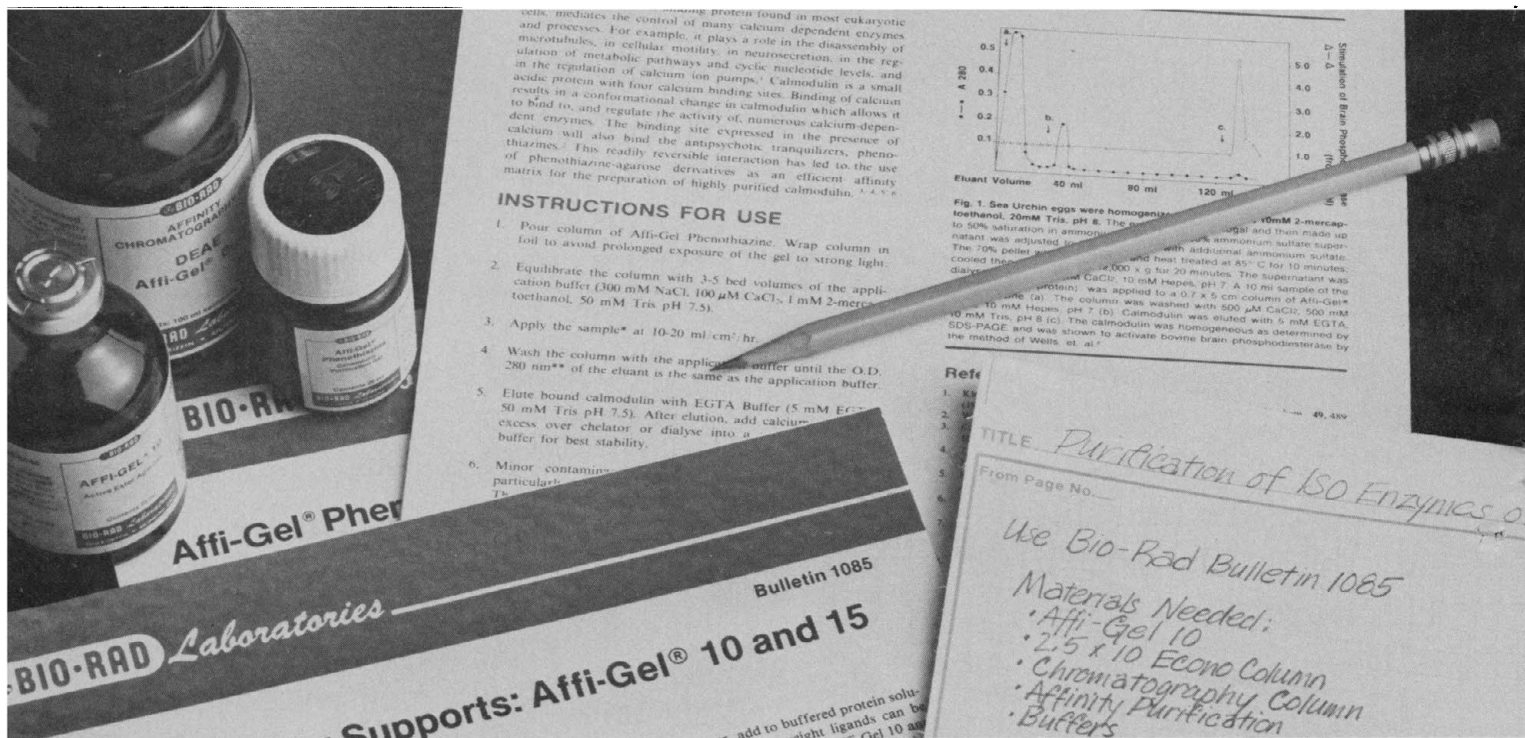

How to stay in step with affinity chromatography

To help you take advantage of today's most powerful and elegant separations techniques, Bio-Rad has assembled an anthology of protocols for the successful use of our Affi-Gel® media for a wide range of applications.

During the past 10 years we have developed practical and reliable step-by-step procedures for each of our Affi-Gel supports.

Now they're compiled in one convenient quick reference folder along with applications data, sample results, and up to date bibliographies. A copy is yours for the asking. We'll also update your anthology as new materials and applications become available.

All Bio-Rad Affi-Gel media come ready to use with a minimum of preparation. They are highly specific materials with excellent stability so you can regenerate and re-use them. Affinity coupling with our Affi-Gel media is more than just a powerful technique—it will save you time and money over conventional purification methods.

For your copy of Bio-Rad's Anthology of Affi-Gel® Protocols, contact:

Circle No. 225 on Readers' Service Card

BIO-RAD

BIO-RAD LABORATORIES
2200 Wright Avenue • Richmond, CA 94804 USA
Phone (415) 234-4130

Also in Rockville Centre, N.Y., Australia, Austria, England, Germany, Italy, Japan, The Netherlands and Switzerland.

AMERICAN ASSOCIATION FOR THE ADVANCEMENT OF SCIENCE

Science serves its readers as a forum for the presentation and discussion of important issues related to the advancement of science, including the presentation of minority or conflicting points of view, rather than by publishing only material on which a consensus has been reached. Accordingly, all articles published in *Science*—including editorials, news and comment, and book reviews—are signed and reflect the individual views of the authors and not official points of view adopted by the AAAS or the institutions with which the authors are affiliated.

Editorial Board

1982: WILLIAM ESTES, CLEMENT L. MARKERT, JOHN R. PIERCE, BRYANT W. ROSSITER, VERA C. RUBIN, MAXINE F. SINGER, PAUL E. WAGGONER, ALEXANDER ZUCKER

1983: FREDERICK R. BLATTNER, BERNARD F. BURKE, CHARLES L. DRAKE, ARTHUR F. FINDEIS, E. PETER GEIDUSCHEK, GLYNN ISAAC, MILTON RUSSELL, WILLIAM P. SLICHTER, JOHN WOOD

Publisher

WILLIAM D. CAREY

Associate Publisher: ROBERT V. ORMES

Editor

PHILIP H. ABELSON

Editorial Staff

Assistant Managing Editor: JOHN E. RINGLE

Production Editor: ELLEN E. MURPHY

Business Manager: HANS NUSSBAUM

News Editor: BARBARA J. CULLITON

News and Comment: COLIN NORMAN (deputy editor), WILLIAM J. BROAD, LUTHER J. CARTER, CONSTANCE HOLDEN, ELIOT MARSHALL, R. JEFFREY SMITH, MARJORIE SUN, JOHN WALSH

Research News: ROGER LEWIN (deputy editor), RICHARD A. KERR, GINA KOLATA, JEAN L. MARX, THOMAS H. MAUGH II, ARTHUR L. ROBINSON, M. MITCHELL WALDROP

Administrative Assistant, News: SCHERRAINE MACK; Editorial Assistants, News: FANNIE GROOM, CASSANDRA WATTS

Senior Editors: ELEANORE BUTZ, MARY DORFMAN, RUTH KULSTAD

Associate Editors: SYLVIA EBERHART, CAITILIN GORDON, LOIS SCHMITT

Assistant Editors: MARTHA COLLINS, STEPHEN KEPPLER, EDITH MEYERS

Book Reviews: KATHERINE LIVINGSTON, Editor; LINDA HEISERMAN, JANET KEGG

Letters: CHRISTINE GILBERT

Copy Editor: ISABELLA BOULDIN

Production: NANCY HARTNAGEL, JOHN BAKER; ROSE LOWERY; HOLLY BISHOP, ELEANOR WARNER; JEAN ROCKWOOD, LEAH RYAN, SHARON RYAN, ROBIN WHYTE

Covers, Reprints, and Permissions: GRAYCE FINGER, Editor; GERALDINE CRUMP, CORRINE HARRIS

Guide to Scientific Instruments: RICHARD G. SOMMER
Assistants to the Editors: SUSAN ELLIOTT, DIANE HOLLAND

Membership Recruitment: GWENDOLYN HUDDLE

Member and Subscription Records: ANN RAGLAND
EDITORIAL CORRESPONDENCE: 1515 Massachusetts Ave., NW, Washington, D.C. 20005. Area code 202. General Editorial Office, 467-4350; Book Reviews, 467-4367; Guide to Scientific Instruments, 467-4480; News and Comment, 467-4430; Reprints and Permissions, 467-4483; Research News, 467-4321. Cable: Advancesci, Washington. For "Information for Contributors," write to the editorial office or see page xi, *Science*, 18 December 1981.

BUSINESS CORRESPONDENCE: Area Code 202. Membership and Subscriptions: 467-4417.

Advertising Representatives

Director: EARL J. SCHERAGO

Production Manager: GINA REILLY

Advertising Sales Manager: RICHARD L. CHARLES

Marketing Manager: HERBERT L. BURKLUND

Sales: NEW YORK, N.Y. 10036: Steve Hamburger, 1515 Broadway (212-730-1050); SCOTCH PLAINS, N.J. 07076: C. Richard Callis, 12 Unami Lane (201-889-4873); CHICAGO, ILL. 60611: Jack Ryan, Room 2107, 919 N. Michigan Ave. (312-337-4973); BEVERLY HILLS, CALIF. 90211: Winn Nance, 111 N. La Cienega Blvd. (213-657-2772); DORSET, VT. 05251: Fred W. Dieffenbach, Kent Hill Rd. (802-867-5581).

ADVERTISING CORRESPONDENCE: Tenth floor, 1515 Broadway, New York, N.Y. 10036. Phone: 212-730-1050.

Essential Federal Information Services

Efficient storage and retrieval of reliable information are of crucial importance to science, medicine, and technology. Of special value has been the development of electronic databases accessible for on-line searches. A decade ago, comprehensive access to scientific information was becoming discouragingly difficult and time-consuming. Now the trend is toward much easier access to the literature. Providing such services is increasingly a function of private enterprise, but two government activities, the National Library of Medicine (NLM) and the National Technical Information Service (NTIS), have had crucial roles which should be continued.

The National Library of Medicine and its predecessors have prepared indexes of the world's medical literature for more than 100 years. In the late 1950's, NLM recognized that the usefulness of its bibliographic publications was being threatened by the increasing volume of literature to be processed. Steps were taken to place the information in machine-readable form. In the early 1960's, NLM was a pioneer in developing equipment for photocomposition. Again, in 1971, NLM pioneered with on-line access to its bibliographic base. Today, NLM indexes more than 20,000 articles a month from about 3000 leading journals. In MEDLINE, access is possible through authors, index terms, or words in titles or abstracts. The databases provided by NLM are highly regarded by research scientists and physicians. They provide an easy way to tap the information produced by the world's biomedical scholars.

Soon after the NLM demonstrated the feasibility of on-line searching of databases in the early 1970's, commercial companies began to develop additional on-line services. The number of records available on-line has expanded sixfold.

The National Technical Information Service facilitates identification and access to specialized information developed by the executive branch of the federal government. Agencies produce reports both on paper and on magnetic tapes. NTIS provides a machine-readable database catalog of these products. Efforts that have led to the reports cost billions of dollars and copies of the reports are sought by many users, especially in the private sector. NTIS has annual sales of about \$25 million and handles more than 25,000 reports and information products per day. Much of the material inventoried by NTIS has limited sales, but a relatively complete file is still maintained. In spite of this, the revenues from sales are such that NTIS is self-supporting.

In keeping with its policy of trying to reduce the size of government, the Reagan Administration has taken steps that threaten the future vitality of the government's two major information services. In spite of the fact that NTIS is self-supporting, disposal of it to private enterprise has been proposed. Private enterprise, unless tightly governed by contract provisions, would be likely to handle only those reports whose sales were profitable.

The NLM has been under budgetary pressure which limits its ability to continue its traditions of innovation. This is unfortunate because the sums of money involved are small in comparison to the service rendered and to the \$200 billion a year or more spent on medical care. NLM has a statutory responsibility to maintain a national medical library and to make its contents effectively available through indexing. Once the index has been prepared, the additional cost of providing on-line access to it is about \$8 million a year, almost all of which is recovered through user charges. Nevertheless, NLM has been the target of an intense lobbying effort by a foreign company aimed at forcing it to increase its charges sharply. This would diminish the value of this great institution to U.S. and world medicine. NLM has achieved an exemplary record in serving its functions. It should be encouraged, not hindered, in its efforts.—PHILIP H. ABELSON

Our freezers are so much in demand because we demand so much of our freezers.

Only a Marine recruit in boot camp goes through more inspections than a Revco laboratory freezer.

After all, when you've been building as many freezers as we have, for as many years as we have, you've got quite a reputation to uphold. Which may explain why our quality control people take their jobs so seriously.



In addition, Revco is the only company in the industry whose *products* — not just a few electrical components — are listed by Underwriters Laboratories, Inc.

We not only establish rigid engineering and design standards, but we also make sure that those standards are applied at every phase of production, in every freezer we make.

And that's the way it should be. Because in work as important as yours, it's imperative for the lab equipment you use to meet high standards. And for more than 40 years, Revco has been the standard by which the industry measures performance and quality.

So send for our new catalog on our complete line of chest and upright freezers. They're available from 3 to 20 cubic feet and from 0°C to -100°C. And we're sure they'll live up to *your* demands, because they've already lived up to *ours*.

REVCO
REVCO
REVCO
REVCO
REVCO
REVCO
REVCO
REVCO

The World's
Leader in
Ultra-Low
Temperature
Equipment

Rheem Manufacturing Company
Scientific Products Division
Aiken Road Rt.1 Box 275
Asheville, N.C. 28804
Telephone 704 658-2711

Call for Symposium Proposals

DETROIT MEETING — 26-31 May 1983

AAAS members are invited to submit symposium proposals for the next Annual Meeting in Detroit, 26-31 May 1983. Please complete the form below, attach a "Synopsis of Objectives" (about 200 words), and send it to us **not later than 1 August 1982**.

We are particularly interested in symposia dealing with the latest developments in science and technology, and the implications of these developments for society.

All symposium proposals are subject to review. If the information submitted is inadequate for reviewing, the proposal will be returned. Endorsement (sponsorship) by a AAAS Section Committee expedites the review process. It is therefore in the interest of the proposer to send a *copy* of the proposal to the appropriate Section Secretary (see table of contents page of

Science for names) for endorsement at the same time the *original* is sent to the AAAS Meetings Office.

Speakers should *not* be confirmed at this time; however, sufficient information about probable speakers and their topics should be provided to allow for evaluation of the proposal.

Some Deadlines

Early October: You will be notified about acceptance, conditional acceptance, or non-acceptance of your proposal. Further information will be provided at that time.

Early November: Preliminary programs with confirmed speakers are due.

Mid January: Final program copy, suitable for publication, is due.

 Annual Meeting Detroit 26 - 31 May 1983	PROPOSED SYMPOSIUM Submit not later than 1 August 1982	Return this form to: AAAS Meetings Office 1101 Vermont Ave., N.W. Washington, D.C. 20005 Telephone: (202) 842-9530
--	---	--

PRINCIPAL ARRANGER (who will receive correspondence):

TELEPHONE _____
(Area Code) (Number)

NAME _____

AFFILIATION _____

MAILING ADDRESS _____
(Street) (City & State) (Zip Code)

CO-ARRANGER (if applicable):

TELEPHONE _____
(Area Code) (Number)

NAME _____

AFFILIATION _____

MAILING ADDRESS _____
(Street) (City & State) (Zip Code)

TITLE OF SYMPOSIUM _____

SPONSORSHIP (AAAS Section, Office, or Affiliate **only**) _____

No. of 3-hr.
Sessions _____

1. Speaker _____

Affiliation _____

Topic _____

2. Speaker _____

Affiliation _____

Topic _____

3. Speaker _____

Affiliation _____

Topic _____

4. Speaker _____

Affiliation _____

Topic _____

5. Speaker _____

Affiliation _____

Topic _____

6. Speaker _____

Affiliation _____

Topic _____

**PLEASE ATTACH SYNOPSIS OF OBJECTIVES OF
PROPOSED SYMPOSIUM (APPROXIMATELY 200 WORDS).**