


SCIENCE

21 September 1956

Volume 124, Number 3221

Editorial	Teacher Certification	513
Articles	Mechanics of Freezing in Living Cells and Tissues: <i>H. T. Meryman</i>	515
	Federal Support of Research in the Life Sciences: <i>W. V. Consolazio and M. C. Green</i>	522
	C. M. Louttit, Psychologist: <i>J. W. Carter, Jr.</i>	526
	H. B. Williams, Physician and Physiologist: <i>W. S. Root, R. H. Kruse, K. S. Cole</i>	527
News of Science	Soviet Physiology; Mathematics Teaching Essay Contest; Bacteriophages in the American Type Culture Collection; Small Animals in the Classroom; U.S.-Brazilian Uranium Agreements; NSF Aids Dissemination of Federal Research Results; Salaries for Chemists; Mauna Loa Observatory; News Briefs; Scientists in the News; Recent Deaths; Education; Miscellaneous	528
	AAAS Building Dedication	534
Reports	Immunologic Comparison of Isolated Surface Membranes of <i>Bacillus megaterium</i> : <i>J. W. Vennes and P. Gerhardt</i>	535
	Algae as Sources of Lysine and Threonine in Supplementing Wheat and Bread Diets: <i>J. M. Hundley, R. B. Ing, R. W. Krauss</i>	536
	San Augustin Plains—Pleistocene Climatic Changes: <i>K. H. Clisby and P. B. Sears</i>	537
	San Augustin Plains—the Geologic Setting: <i>C. E. Stearns</i>	539
	San Augustin Plains—the Sediments: <i>F. Foreman</i>	539
Book Reviews	<i>Noradrenaline; The Changing Universe; New Lives for Old; The Pursuit of Science in Revolutionary America, 1735–1789; Techniques Générales du Laboratoire de Physique</i> ; New Books; Miscellaneous Publications	540
Meetings and Societies	New York Meetings of the AAAS, 1887–1956; Meeting Notes; Forthcoming Events	544

QUALITY, DEPENDABILITY, PERFORMANCE

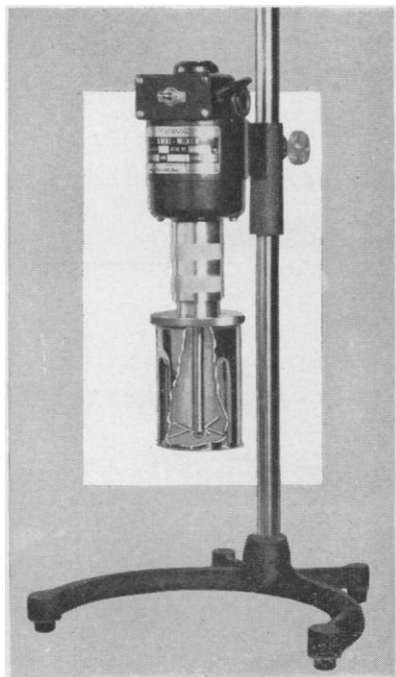


OMNI - MIXER: A portable, high-speed instrument designed exclusively by Ivan Sorvall, Inc. for all laboratory homogenizing, mixing, emulsifying and difficult disintegrating work. Specially-shaped, stainless-steel, 200-ml mixing chamber easily cooled or heated. Accommodates mason jars • rugged construction • easily cleaned and sterilized.

Write for Bulletin SC-96

ULTRA-MICROTOME: 1/40 to 1/2 Micron. Simple to operate—makes Serial Sectioning for Electron Microscopy a routine laboratory technique. Mechanical advancement • collet-type chuck and many other features • Lightweight and portable.

Write for Bulletin SC-95

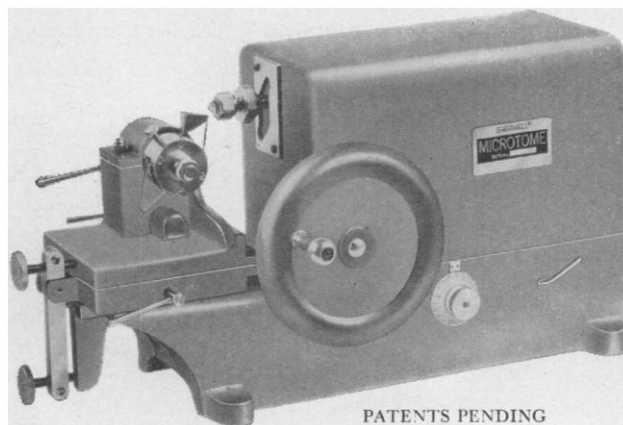

PATENTS PENDING



SUPERSPEED ANGLE CENTRIFUGE: 15,500-rpm—31,000 x G—400-ml. Portable, with separate transformer for remote control—may be operated in cold room, incubator, etc. Accommodates tubes from 1- to 50-ml (with use of adapters). Self-balancing, direct drive for smooth, trouble-free operation. "The workhorse of the modern research laboratory."

Write for Bulletin SC-91


M, SP, & SP/X ANGLE CENTRIFUGES: 5,000-rpm—3,500 x G.
Type M—Twenty-Seven 15-ml Tubes
Type SP—Six 75-ml Tubes
Type SP/X—Ten 50-ml Tubes and Five 15-ml Tubes. Tube adapters available for greater versatility. Other sizes of Angle Centrifuges available.

Write for Bulletin SC-91



SUPERSPEED REFRIGERATED CENTRIFUGE: 15,500-rpm—31,000 x G. Exclusive dual automatic temperature and holding control establishes and maintains temperatures down to and below 0°C. No external pre-cooling of rotors is necessary. A large selection of "Servall" rotors, in addition to the Superspeed Rotor, may be interchanged easily and rapidly on the "Servall"-originated direct drive. No complicated drive attachments are required in order to obtain the higher speeds. Maximum capacity at Superspeeds: 400-ml Maximum capacity at lower speeds: 2,000-ml The "Servall" is the least expensive Refrigerated Centrifuge in the Superspeed Range.

Write for Bulletin SC-91R

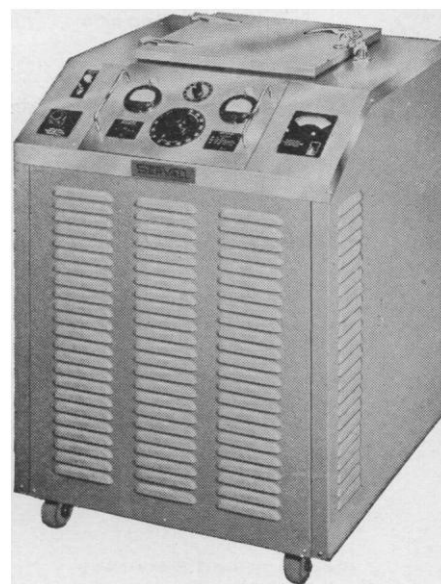

ALSO, IVAN SORVALL, INC. ARE SOLE DISTRIBUTORS IN THE UNITED STATES AND CANADA FOR LKB PHYSICAL CHEMISTRY INSTRUMENTS.

Ivan Sorvall, Inc.

NORWALK • CONNECTICUT

Designers, Manufacturers and Distributors of "Servall" Specialized Laboratory Instruments Since 1930.

Rapid centrifugal preparation

with the *Spinco* Model L
ULTRACENTRIFUGE

- **VACUUM**
- **162 ml – 40,000 rpm**
- **REFRIGERATED**



FULLY-AUTOMATIC features make the Spinco Model L Preparative Ultracentrifuge operationally self-supervising. Simple speed- and time-controls are preset by operator and electronic unit brings rotor to speed, times run, and decelerates with smoothly operating automatic brake.

Various sizes and types of interchangeable rotors mount in operating chamber with instantaneous coupling. Time between runs, using pre-loaded rotors, is less than a minute. Speed and vacuum are established almost simultaneously. Acceleration time for the 40,000-rpm rotor is only four minutes.

SAFETY FEATURES include automatic overspeed cutout; interlocks to prevent running at high speeds without vacuum in chamber and prevent breaking vacuum or opening chamber while running at high speeds. Rotor chamber is encased in thick double steel guard.

UNIT IS SELF-CONTAINED, of worktable height, and requires only electric connection for installation. Spinco electric drive and gearing system are of types proved in more than 100,000 hours of actual field operation. Self-balancing feature eliminates weighing of loaded tubes. Controllable refrigeration holds rotors and contents at any temperature down to -5°C .

In operation, the Model L applies forces up to 144,000 times gravity on 162 ml of material at a maximum speed of 40,000 rpm. Individual rotor tubes have sealing caps—are available in plastic or metal. Rotors provide second vacuum seal in quickly-detachable cover with integral lifting device. Wide selection of rotors is available with largest holding 1700 ml—

*Send for
details*

Spinco *division*

BECKMAN INSTRUMENTS, INC.
BELMONT 5, CALIFORNIA

AAAS Building Dedication

To celebrate the completion of its new headquarters building, the Association will hold a special meeting in Washington, D.C., on 12 October 1956. In keeping with the nature of the Association, the program will be devoted to a scientific topic: "The Uses and Effects of Atomic Radiation." All members of the Association and of the interested public are cordially invited. The meeting will be held at the Carnegie Institution of Washington, 16th and P Streets, NW.

Program

THE USES AND EFFECTS OF ATOMIC RADIATION

Chairman: Paul B. Sears, President

10 A.M.

Welcome

Caryl P. Haskins, *President, Carnegie Institution of Washington*

Radiation and the Human Body

Shields Warren, *Pathologist, New England Deaconess Hospital*

Radiation and Genetics

L. C. Dunn, *Professor of Zoology, Columbia University*

The Uses of Atomic Radiation and Energy

Lawrence R. Hafstad, *Vice President, General Motors Corporation*

2 P.M.

Social Implications of Atomic Radiation and Energy

Detlev W. Bronk, *President, Rockefeller Institute*

Current Research Findings

Willard F. Libby, *Commissioner, U.S. Atomic Energy Commission*

What We Most Need To Know

Laurence H. Snyder, *President-Elect of the Association and
Dean of the Graduate College, University of Oklahoma*

Dedication of the American Association for the Advancement of Science Building

Paul B. Sears, *President of the Association and
Chairman, Conservation Program, Yale University*

At the conclusion of the meeting, the Association will hold open house at its new building at 1515 Massachusetts Avenue, NW.