

# SCIENCE

28 September 1956

Volume 124, Number 3222


<b>Editorial</b>	Manpower and Education .....	565
<b>Articles</b>	Scientists in the British Civil Service: <i>E. McCrensky</i> .....	567
	Uptake and Turnover of Calcium-45 by the Guppy: <i>H. L. Rosenthal</i> .....	571
	A. H. Hersh, Geneticist: <i>D. P. Quiring</i> and <i>F. DeMarinis</i> .....	575
<b>News of Science</b>	Indians in the Peruvian Andes; How Many Chromosomes? International Co-operation in Uranium Exploration; National Seed Storage Laboratory; U.N. Exhibit to Commemorate Geneva Nuclear Conference; Archeology in Alabama; Organized Labor and the New Michigan Reactor; Borneo Zoological Expedition; TB in the United States; British Research Budget; Abortion in the U.S.S.R.; Radiation in Monkeys; News Briefs; Scientists in the News; Recent Deaths; Education; Grants, Fellowships, and Awards; In the Laboratories; Miscellaneous .....	576
<b>Reports</b>	Subjective Evaluation and Reinforcing Effect of a Verbal Stimulus: <i>G. Mandler</i> and <i>W. K. Kaplan</i> .....	582
	Influence of "Aging" on the Characteristics of an Electrodeless Discharge: <i>C. N. R. Rao</i> .....	583
	Preparation of Pure Hydroxyapatite Crystals: <i>A. Perloff</i> and <i>A. S. Posner</i> ....	583
	Beryllium-10 Produced by Cosmic Rays: <i>J. R. Arnold</i> .....	584
	Quantum Yields of Fluorescence of Plant Pigments: <i>P. Latimer</i> , <i>T. T. Bannister</i> , <i>E. Rabinowitch</i> .....	585
	Sperm Transport in the Reproductive Tract of the Female Rabbit: <i>G. S. Greenwald</i> .....	586
	Cold-Adapted Genetic Variants of Polio Viruses: <i>G. R. Dubes</i> and <i>M. Chapin</i> ..	586
	Influence of Crystalline Elastase on Experimental Atherosclerosis in the Chicken: <i>D. M. Tennent</i> et al. ....	588
<b>Letters</b>	Scientific Method, Statistical Inference, and the Law: <i>G. W. McElrath</i> and <i>J. E. Bearman</i> .....	589
	Of Books and Reading: <i>L. Wallerstein</i> ; <i>G. Gallup</i> .....	590
<b>Book Reviews</b>	<i>Reactor Shielding Design Manual</i> ; <i>Floods</i> ; <i>Methods of Biochemical Analysis</i> ; <i>Studies of the Psychology and Behavior of Captive Animals in Zoos and Circuses</i> ; <i>Précis de Géologie</i> ; <i>The Negro Potential</i> ; Books Reviewed in <i>The Scientific Monthly</i> ; New Books .....	591
<b>Meetings and Societies</b>	Election of AAAS Officers; Physiology of Reproduction; Meeting Notes; Society Elections; Forthcoming Events .....	594

SARGENT-MALMSTADT

# Automatic Titrator

Using Second Derivative Method of H. V. Malmstadt

1. Greater operational time saving than any former so-called automatic equipment.
2. No electrode problems. Most titrations can be carried out with plain platinum or platinum-rhodium wire electrodes.
3. Correct end point is automatically derived by the circuit — no statistical or theoretical studies needed. No pre-setting of volt-
- ages. Push-button start, automatic shutoff.
4. Inherent accuracy is superior to voltage measurement or color indication.
5. Precision better than limits of normal techniques and reading.
6. Directly applicable to all types of redox, acid-base, precipitation and complexation titration systems.
7. Low cost, minimum maintenance.

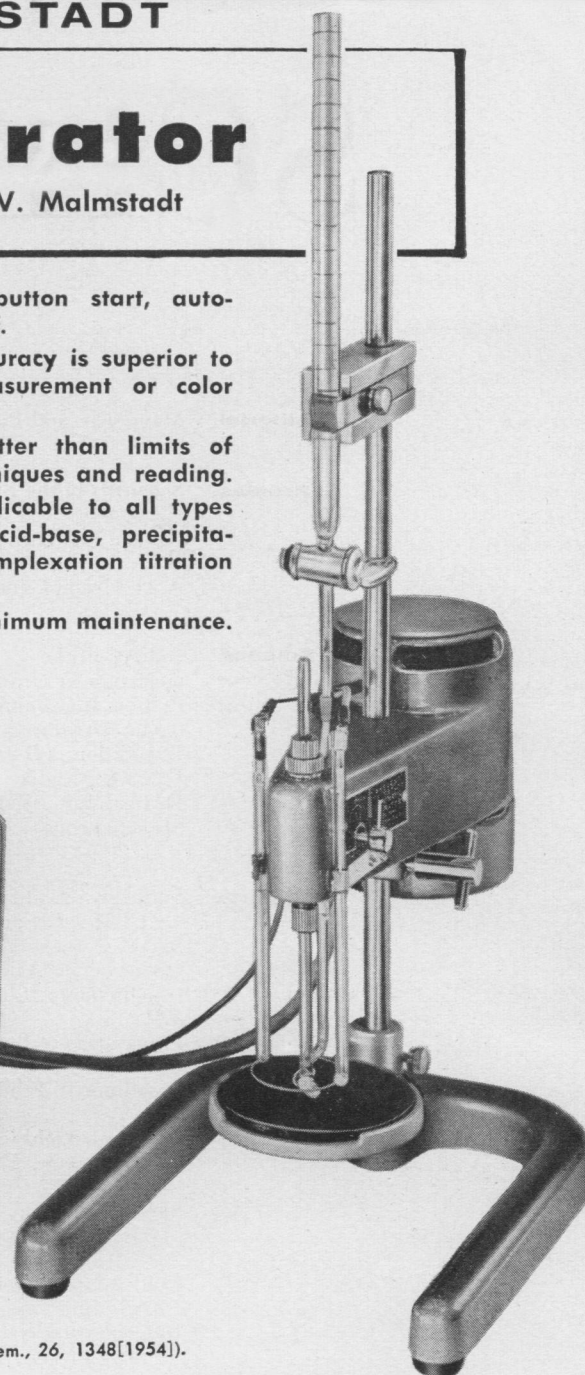

Developed and manufactured by E. H. Sargent & Co.,  
based on the original work of H. V. Malmstadt, (Anal. Chem., 26, 1348[1954]).

**S-29690 ELECTROMETRIC TITRATION APPARATUS**—Automatic, Potentiometer, Sargent-Malmstadt (Patents pending).

For use in all titrations usually classified as potentiometric. Consisting fundamentally of a mixing and delivery unit with associated supports and an electronic control circuit.

The control unit contains the detecting, amplifying, differentiating and relaying circuits with a simple operational panel. At the rear of the control unit are a stirring motor switch for manual use, a fuse, connecting outlets to the valve unit, connecting outlet to the power supply, and a recorder outlet permitting, for research purposes, the recording of the second derivative curve with automatic cutoff.

The delivery unit contains a valve mechanism with an adjustable rate controller, an open ended motor driven chuck for the stirrer, stirring motor and electrode holders. The valve mechanism,

mounted on a support with swinging beaker platform and a clamp for the titrating burette, is designed to accommodate any form of burette with or without stopcocks, so that it may apply to any type of titrating installation.

Dimensions of control unit: height, 7 $\frac{3}{4}$  inches; width, 10 inches; depth, 9 inches; power consumption, 100 watts; total shipping weight, 35 lbs.

Complete, including: control unit; stirring valve unit; support stand with beaker platform and burette clamp; one S-30490 calomel electrode; one each S-30440 electrode, platinum and platinum rhodium; one S-29695 delivery tip with internal platinum electrode, one S-76667 stirring rod; one S-76668 large stirring rod; connecting cable for recorder; and connecting cord and plug for standard outlets. Without beakers or burettes. For operation from 115 vol, 60 cycle A.C. circuits..... **\$450.00**

## SARGENT

SCIENTIFIC LABORATORY INSTRUMENTS • APPARATUS • SUPPLIES • CHEMICALS

*Catalog No 100  
Now Available*

E. H. SARGENT & COMPANY, 4647 W. FOSTER AVE., CHICAGO 30, ILLINOIS  
MICHIGAN DIVISION, 8560 WEST CHICAGO AVENUE, DETROIT 4, MICHIGAN  
SOUTHWESTERN DIVISION, 5915 PEELER STREET, DALLAS 19, TEXAS  
SOUTHEASTERN DIVISION, 3125 SEVENTH AVE., N., BIRMINGHAM 4, ALA.

# NOW

# E-X-T-E-N-D-E-D

# Range!



## NEW TESTS, NEW APPLICATIONS

Simple, inexpensive tests, helpful in diagnosis of coronary artery disease, can now be made with Bausch & Lomb Spectronic 20 Colorimeters. Special adaptation, extending instrument range to  $340m\mu$  in the U-V, makes assay of serum glutamic-oxalacetic transaminase a routine clinical procedure. (Detailed clinical method available on request.)

Many other new tests, including vitamin assays, are made possible by pinpoint selection of any wavelength in the  $340m\mu$ - $950m\mu$  range.

## SAME EASE, SPEED AND DEPENDABLE ACCURACY


*Easy As Tuning Your Radio!* Dial instantly sets Certified-Precision diffraction grating to desired wavelength; no color filters to fuss with.

*Fastest Readings!* Simply insert sample tube; instant-acting meter gives exact percent transmission, or optical density.

*Accurate Analyses!* Narrow band pass (only  $20m\mu$ ) assures highest spectral purity. 49 clinical calibrations available.

## LOW PRICE, DOUBLE VALUE

Colorimeter *plus* spectrophotometer,  $375m\mu$ - $950m\mu$  (extended range to  $340m\mu$  at slight extra cost), for less than the price of most colorimeters alone!

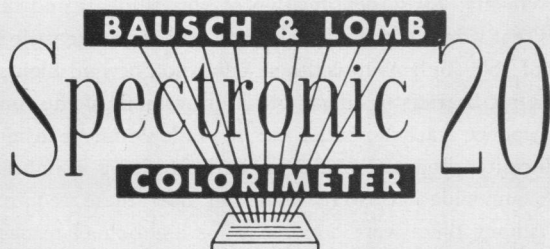

**BAUSCH & LOMB OPTICAL CO.**  
64209 ST. PAUL ST., ROCHESTER 2, N. Y.

- ☐ Send me 11-page reprint of "Clinical Method for Assay of Serum Glutamic-Oxalacetic Transaminase" plus references and charts.
- ☐ Send me B&L Spectronic 20 Catalog D-266.
- ☐ I would like an obligation-free demonstration of the B&L Spectronic 20 Colorimeter at my convenience.

NAME .....

TITLE .....

ADDRESS .....

CITY ..... ZONE ..... STATE.....

**MAIL COUPON TODAY  
FOR VALUABLE DATA**

**BAUSCH & LOMB**



America's only complete optical source... from glass to finished product.